

March 2020

artus[®] JCV RG PCR Kit

Instructions for Use (Handbook)



Version 1

For use with Rotor-Gene[®] Q instruments

IVD

CE

REF

4532265



altona Diagnostics GmbH,
Mörkenstraße 12, 22767 Hamburg, GERMANY

R3 MAT

1121142-EN

Distributed by QIAGEN GmbH, QIAGEN Strasse 1, 40724 Hilden, GERMANY

Contents

Intended Use	4
Summary and Explanation	4
Pathogen information	4
Principle of the Procedure	5
Materials Provided	6
Kit contents	6
Materials Required but Not Provided	7
Warnings and Precautions	8
Warnings	8
Precautions	8
Reagent Storage and Handling	10
Kit components	10
Procedure	12
DNA extraction	12
Protocol: Detection of JCV-specific DNA	14
Interpretation of Results	24
Run validity	24
Qualitative analysis	25
Quantitative analysis	26
Limitations	28
Quality Control	28
Performance Characteristics	29

Analytical sensitivity.....	29
Analytical specificity	30
Linear range	31
Precision	31
Repeatability.....	32
Symbols.....	34
Troubleshooting Guide	36
Ordering Information	37
Handbook Revision History	41

Intended Use

The *artus*® JCV RG PCR Kit (96) is an *in vitro* diagnostic test, based on real-time PCR technology, for the detection and quantification of JC Virus (JCV) specific DNA.

Summary and Explanation

The *artus* JCV RG PCR Kit constitutes a ready-to-use system for the detection of JCV-specific DNA using real-time PCR on Rotor-Gene Q instruments. The assay includes a heterologous amplification system (Internal Control) to identify possible PCR inhibition and to confirm the integrity of the kit reagents.

Pathogen information

The family *Polyomaviridae* consists of at least 16 members, infecting different mammalian species. Two of them, the human polyomaviruses BK virus (BKV) and JC virus (JCV), establish ubiquitous infections worldwide. Primary infections generally occur in early childhood and are typically subclinical, followed by a lifelong persistence. The epithelial cells of kidney, lymphocytes, and oligodendrocytes were identified as predominant cell types persistently infected by JCV.

Whereas polyomavirus infections remain unapparent in immunocompetent individuals, reactivation of BKV and/or JCV in association with immunosuppression may lead to serious diseases. JCV reactivation may result in development of progressive multifocal leukoencephalopathy (PML), a demyelinating disease of the central nervous system.

Polyomaviruses are small non-enveloped, icosahedral viruses with a supercoiled double-stranded DNA genome of approximately 5000 base pairs. Human polyomaviruses share ~70% sequence homology with simian virus 40 (SV40). Despite this high sequence homology, polyomaviruses exhibit a restricted host range with distinct biological behavior and disease pathogenesis.

Principle of the Procedure

The JCV RG Master A and JCV RG Master B contain reagents and enzymes for the specific amplification of target regions within the JCV genome and for the direct detection of the specific amplicon in the fluorescence channel Cycling Green of Rotor-Gene Q instruments.

In addition, the *artus* JCV RG PCR Kit contains a heterologous amplification system to identify potential failures during the assay process. This is detected as an Internal Control (IC) in fluorescence channel Cycling Yellow of Rotor-Gene Q instruments.

Probes specific for JCV DNA are labeled with the fluorophore FAM™. The probe specific for the Internal Control (IC) is labeled with the fluorophore JOE™. The use of probes labeled with spectrally distinguishable fluorophores enables simultaneous detection and quantification of JCV DNA as well as detection of the Internal Control in the corresponding channels of the Rotor-Gene Q instrument.

Materials Provided

Kit contents

artus JCV RG PCR Kit		(96)
Catalog number		4532265
Number of reactions		96
Blue	JCV RG Master A	8 x 60 µl
Purple	JCV RG Master B	8 x 180 µl
Green	JCV RG IC	1 x 1000 µl
Red	JCV RG QS*	4 x 250 µl
White	H ₂ O	1 x 500 µl
	Handbook	1

*The *artus* JCV RG PCR Kit contains 4 Quantification Standards (QS1–QS4).

Materials Required but Not Provided

Prior to use, ensure that instruments have been checked and calibrated according to the manufacturer's recommendations.

Reagents

- QIAamp DNA Mini Kit (QIAGEN cat. no. 51304 or 51306; See "DNA extraction", page 12)

Consumables

- 0.1 ml Strip Tubes and Caps, for use with 72-well rotor (QIAGEN, cat. no. 981103 or 981106)
- Nuclease-free, low DNA-binding microcentrifuge tubes for preparing master mixes
- Nuclease-free pipet tips with aerosol barriers

Equipment

- Rotor-Gene Q MDx 5plex, Rotor-Gene Q 5plex or Rotor-Gene Q 6plex instrument
- Rotor-Gene Q software version 2.3.1 or higher
- Loading Block 72 x 0.1 ml Tubes, aluminum block for manual reaction setup (QIAGEN, cat. no. 9018901)
- Dedicated adjustable pipets for sample preparation
- Dedicated adjustable pipets for PCR master mix preparation
- Dedicated adjustable pipets for dispensing template DNA
- Vortex mixer
- Benchtop centrifuge with rotor for 2 ml reaction tubes

Warnings and Precautions

For in vitro diagnostic use.

Read all instructions carefully before using the test.

Warnings

When working with chemicals, always wear a suitable lab coat, disposable gloves, and protective goggles.

Precautions

- Use of this product is limited to personnel specially instructed and trained in the techniques of real-time PCR and in vitro diagnostic procedures.
- Specimens should always be treated as infectious and/or biohazardous in accordance with safe laboratory procedures.
- Wear protective disposable powder-free gloves, a laboratory coat, and eye protection when handling specimens.
- Avoid microbial and nuclease (DNase/RNase) contamination of the specimen and the components of the kit.
- Always use DNase-/RNase-free disposable pipet tips with aerosol barriers.
- Always wear protective disposable powder-free gloves when handling kit components.
- Use separated and segregated working areas for specimen preparation, reaction setup and amplification/detection activities. The workflow in the laboratory should proceed in a unidirectional manner. Always wear disposable gloves in each area, and change them before entering different areas.
- Dedicate supplies and equipment to the separate working areas, and do not move them from one area to another.

-
- Store positive and/or potentially positive material separately from all other components of the kit.
 - Do not open the reaction tubes/plates post amplification to avoid contamination with amplicons.
 - Additional controls may be tested according to guidelines or requirements of local, state, and/or federal regulations or accrediting organizations.
 - Do not use components of the kit that have passed their expiration date.
 - Discard sample and assay waste according to your local safety regulations.

Reagent Storage and Handling

Kit components

The *artus* JCV RG PCR Kit is shipped on dry ice. The components of the kit should arrive frozen. If one or more components are not frozen upon receipt, or if tubes have been compromised during shipment, contact QIAGEN Technical Services for assistance. Upon receipt, store all kit components at -30°C to -15°C .

Avoid thawing and freezing Master reagents more than two times as this may reduce assay performance. Freeze the reagents in aliquots if they are to be used intermittently. Do not store reagents at 4°C for longer than 2 hours. Protect JCV RG Master A and JCV RG Master B from light.

The *artus* JCV RG PCR Kit includes:

- Two Master reagents (JCV RG Master A and JCV RG Master B)
- Template Internal Control (JCV RG IC)
- Four Quantification Standards (JCV RG QS1–QS4)
- PCR-grade water (H_2O)

JCV RG Master A and JCV RG Master B reagents contain all components (buffer, enzymes, primers, and probes) for the amplification and detection of JCV-specific DNA and the Internal Control in a single reaction.

The Quantification Standards contain standardized concentrations of JCV-specific DNA. These Quantification Standards were calibrated against the 1st WHO International Standard for JC virus (JCV) for Nucleic Acid Amplification Techniques (NIBSC code 14/114). These can be used individually as positive controls or together to generate a standard curve, which can be used to determine the concentration of JCV-specific DNA in the sample. The concentrations of the Quantification Standards are shown in Table 1.

Table 1. Concentration of Quantification Standards

Quantification Standard	Concentration (IU/μl)
QS1	10,000
QS2	1000
QS3	100
QS4	10

Procedure

DNA extraction

JCV-specific target sequences are amplified from DNA. As assay performance is dependent on the quality of the template DNA, make sure to use a sample preparation kit that yields DNA suitable for use in downstream PCR.

The QIAamp DNA Mini Kit (QIAGEN, cat. no. 51304 or 51306) is recommended for DNA purification for use with the *artus* JCV RG PCR Kit. Carry out DNA purification according to the instructions in the *QIAamp DNA Mini Handbook*.

As the wash buffers in the QIAamp DNA Mini Kit contain ethanol, carry out an additional centrifugation step prior to elution. Place the QIAamp Mini spin column in a new 2 ml collection tube and discard the old collection tube with the filtrate. Centrifuge for 10 minutes at approximately 17,000 x *g* (~13,000 rpm) in a benchtop centrifuge.

Important: The use of carrier RNA is crucial for extraction efficiency and stability of the extracted nucleic acid.

Important: Ethanol is a strong inhibitor in real-time PCR. If your sample preparation kit uses wash buffers containing ethanol, make sure to remove all traces of ethanol prior to elution of the nucleic acid.

Internal Control

The *artus* JCV RG PCR Kit contains a heterologous Internal Control, which can either be used as a PCR inhibition control, or as a control of the sample preparation procedure (nucleic acid extraction) and as a PCR inhibition control (step 2a, page 14).

If the Internal Control is used as a PCR inhibition control but not as a control for the sample preparation procedure, add the Internal Control directly to the mixture of JCV RG Master A and JCV RG Master B, as described in step 2b of the protocol (page 15).

Regardless of which method/system is used for nucleic acid extraction, the Internal Control must not be added directly to the specimen. The Internal Control should always be added to the specimen/lysis buffer mixture. The volume of Internal Control to be added to the specimen/lysis buffer mixture depends only on the elution volume and represents 10% of the elution volume. For example, when using the QIAamp DNA Mini Kit, the DNA is eluted in 60 µl Buffer AE. Therefore, add 6 µl Internal Control to the specimen/lysis buffer mixture of each sample.

Important: Do not add the Internal Control and/or the carrier RNA directly to the specimen.

Protocol: Detection of JCV-specific DNA

Important points before starting

- Before beginning the procedure, read “Precautions”, page 8.
- Take time to familiarize yourself with the Rotor-Gene Q instrument before starting the protocol. See the instrument user manual.
- Make sure that at least one positive control and one negative control (PCR-grade water) are included per PCR run.

Things to do before starting

- Make sure that the cooling block (accessory of the Rotor-Gene Q instrument) is pre-cooled to 2–8°C.
- Before each use, all reagents need to be thawed completely, mixed (by repeated up and down pipetting or by quick vortexing) and centrifuged briefly.

Procedure

1. Place the desired number of PCR tubes into the adapters of the cooling block.
2. If you are using the Internal Control to monitor the DNA isolation procedure and to check for possible PCR inhibition, follow step 2a. If you are using the Internal Control exclusively to check for PCR inhibition, follow step 2b.

If the Internal Control was added during the sample preparation procedure, then at minimum, the negative control – which is not a negative specimen – must include the Internal Control.

- 2a. The Internal Control has already been added to the isolation (see “Internal Control”, page 12). In this case, prepare a master mix according to Table 2.

The reaction mix typically contains all of the components needed for PCR, except the sample.

Table 2. Preparation of master mix (Internal Control used to monitor DNA isolation and to check for PCR inhibition)

Component	1 reaction	12 reactions
JCV RG Master A	5 µl	60 µl
JCV RG Master B	15 µl	180 µl
Total volume	20 µl	240 µl

- 2b. The Internal Control must be added directly to the mixture of JCV RG Master A and JCV RG Master B. In this case, prepare a master mix according to Table 3.

The reaction mix typically contains all of the components needed for PCR, except the sample.

Table 3. Preparation of master mix (Internal Control used exclusively to check for PCR inhibition)

Component	1 reaction	12 reactions
JCV RG Master A	5 µl	60 µl
JCV RG Master B	15 µl	180 µl
JCV RG IC	1 µl	12 µl
Total volume	21 µl	252 µl

* The volume increase caused by adding the Internal Control is neglected when preparing the PCR assay. The sensitivity of the detection system is not impaired.

3. Pipet 20 µl of the master mix into each PCR tube. Then add 10 µl of the eluted sample DNA and mix well by pipetting repeatedly up and down. Correspondingly, add 10 µl of a positive control or Quantification Standard or 10 µl H₂O (PCR-grade water) as a negative control.

Make sure to have at least one positive control and one negative control per run. For quantification, use all 4 Quantification Standards (QS1–QS4).

4. Close the PCR tubes. Make sure that the locking ring (accessory of the Rotor-Gene instrument) is placed on top of the rotor.

5. For the detection of JCV-specific DNA, create a temperature profile according to the following steps.

Setting the general assay parameters	Figures 1, 2, 3, 4
Initial activation of the hot-start enzyme	Figure 5
Amplification of the DNA	Figure 6
Adjusting the fluorescence channel sensitivity	Figure 7
Starting the run	Figure 8

All specifications refer to the Rotor-Gene Q software version 2.3.1 and higher. Please find further information on programming Rotor-Gene instruments in the instrument user manual. In the illustrations, these settings are framed in bold black.

6. First, open the **New Run Wizard** dialog box with the **Advanced** version and select **Two Step** (Figure 1). Click **New** to continue.

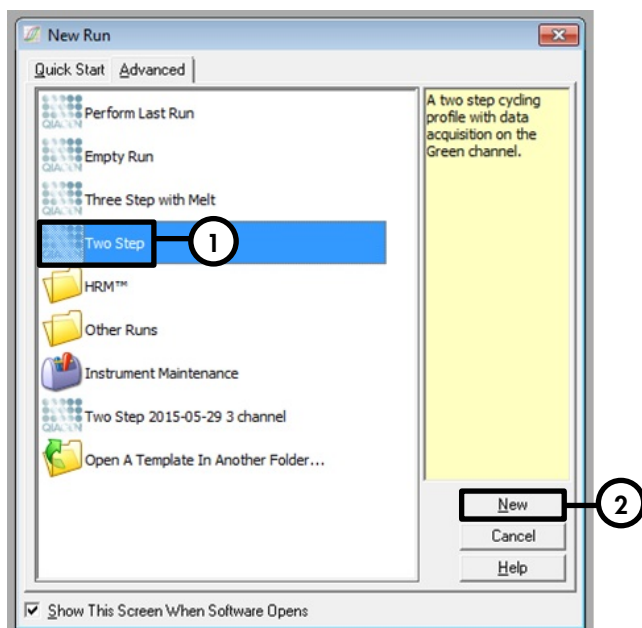


Figure 1. The New Run dialog box.

7. In the next **New Run Wizard** dialog box (Figure 2), check the **Locking Ring Attached** box and click **Next**.

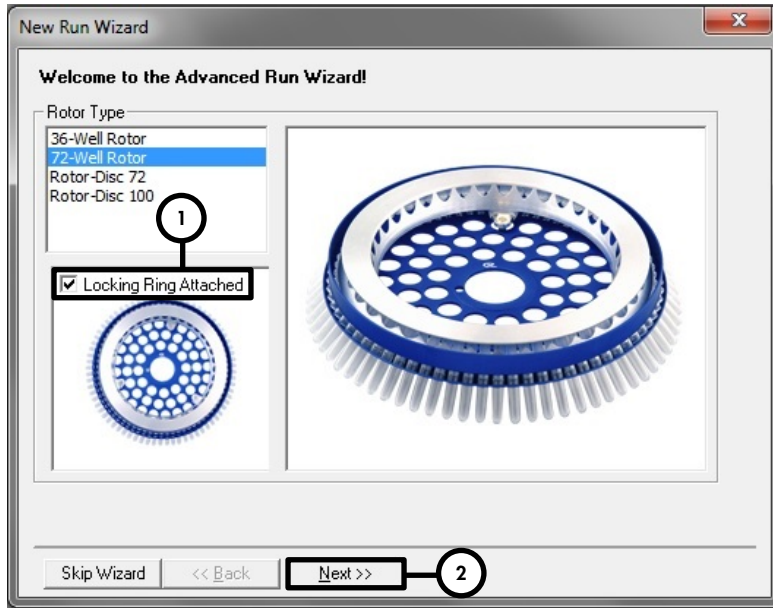


Figure 2. The New Run Wizard dialog box.

8. Select **30** for the PCR reaction volume and click **Next** (Figure 3).

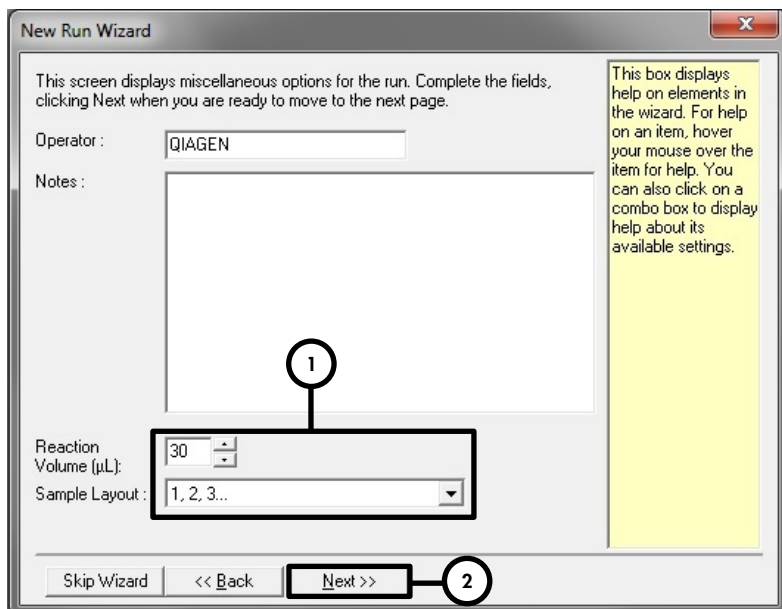


Figure 3. Setting the general assay parameters.

9. Click the **Edit Profile** button in the next **New Run Wizard** dialog box (Figure 4), and program the temperature profile as shown in Figures 5–6.

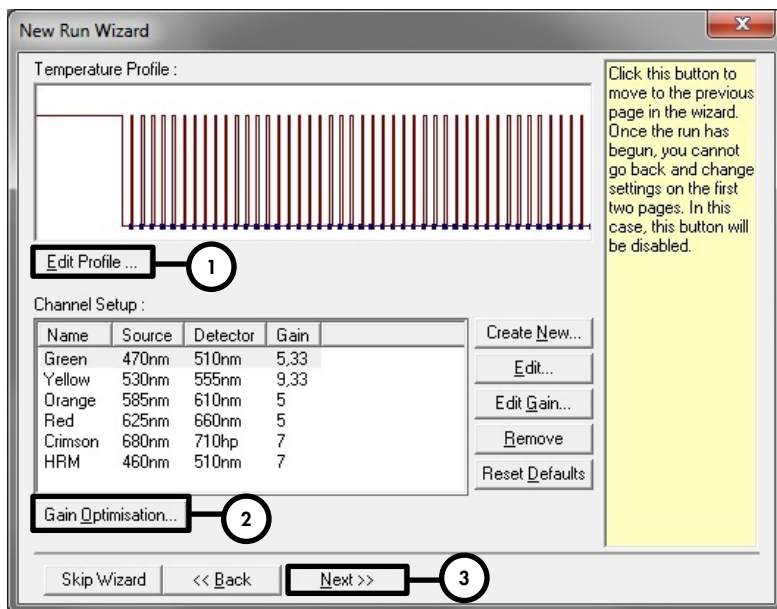
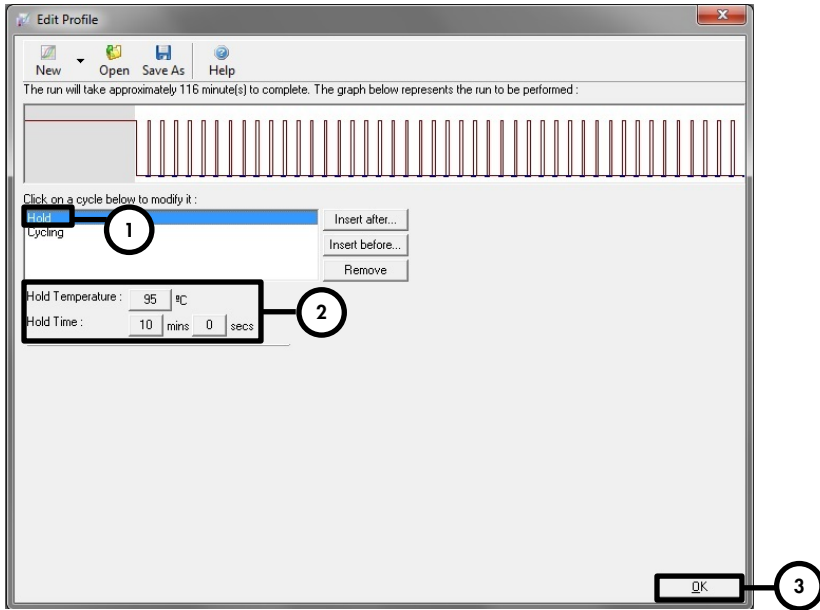


Figure 4. Editing the profile.



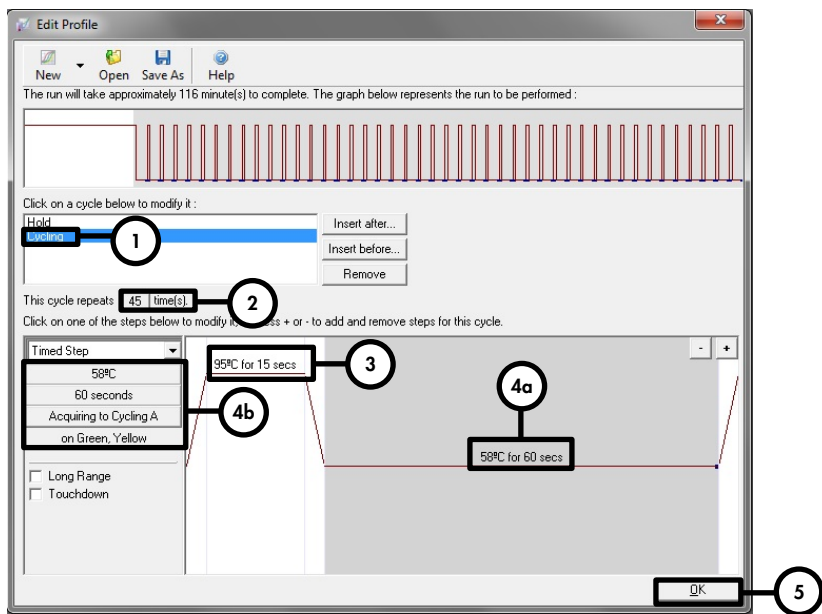


Figure 6. Amplification of the DNA.

10. The detection range of the fluorescence channels has to be determined according to the fluorescence intensities in the PCR tubes. Click **Gain Optimisation** in the **New Run Wizard** dialog box (see Figure 4, Step 2) to open the **Auto-Gain Optimisation Setup** dialog box (Figure 7). Check the **Perform Optimisation Before 1st Acquisition Box** (Figure 7). Make sure that both channels (Green and Yellow) are selected for **Auto-Gain Optimisation** (Figure 7). (Find channels in the drop-down menu under **Channel Settings** and click **Add**.) Click **Close** of the **Auto-Gain Optimisation Setup** dialog box when the gain calibration is completed.

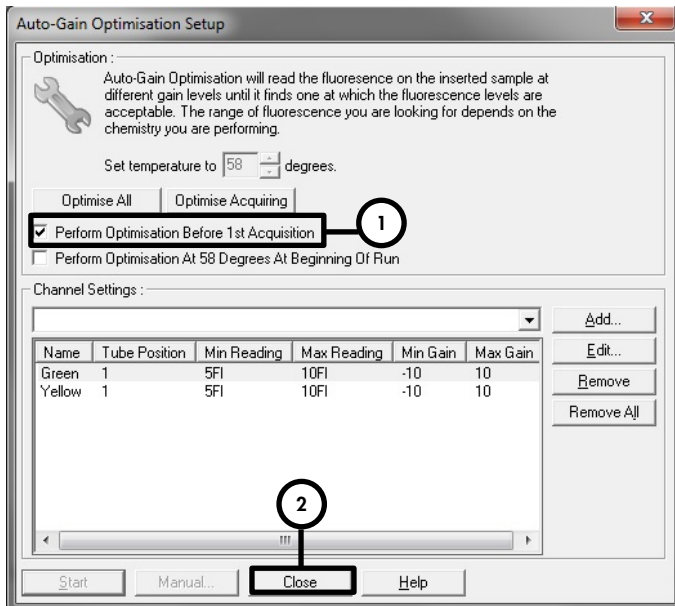


Figure 7. Adjusting the fluorescence channel sensitivity.

11. The gain values determined by the channel calibration are saved automatically and are listed in the last menu window of the programming procedure (Figure 8). Click **Start Run**.

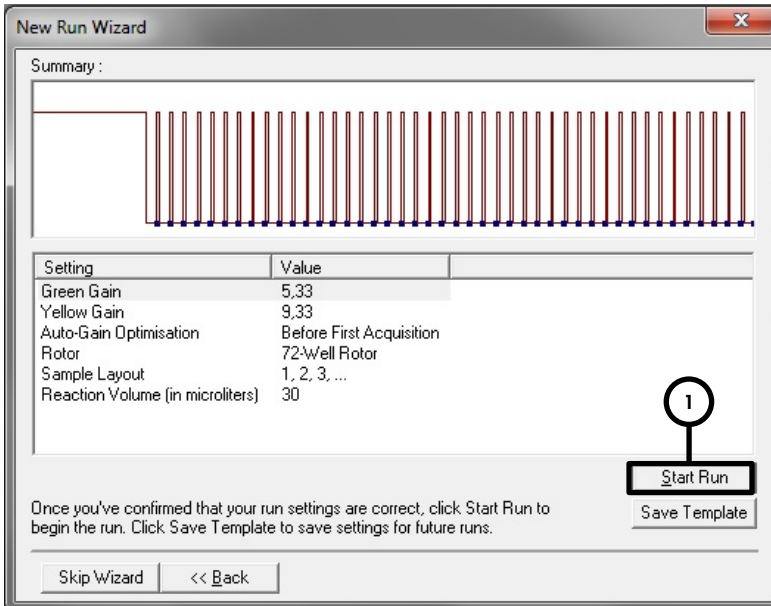


Figure 8. Starting the run.

12. After the run is finished, analyze the data (see "Interpretation of Results", page 24).

Interpretation of Results

Run validity

Valid qualitative run

The following control conditions must be met for a qualitative run to be valid (Table 4).

Table 4. Control conditions for a valid qualitative run

Control ID	Detection channel	
	Cycling Green	Cycling Yellow
Positive control (QS)	POSITIVE	POSITIVE
Negative control	NEGATIVE	POSITIVE

Invalid qualitative run

A qualitative run is invalid if the run has not been completed or if any of the control conditions for a valid qualitative run have not been met.

In case of an invalid qualitative run, repeat the PCR or extract DNA from the original samples again if no DNA is left over.

Valid quantitative run

A quantitative run is valid if all control conditions for a valid qualitative run have been met (see Table 4, above). Furthermore, for accurate quantification results, a valid standard curve needs to be generated. For a valid quantitative run, the standard curve must have the following control parameter values (Table 5).

Table 5. Control parameters for a valid standard curve

Control parameter	Valid value
Slope	-3.743/-2.765
PCR efficiency	85%/130%
R squared (R ²)	>0.98

Invalid quantitative run

A quantitative run is invalid if the run has not been completed or if any of the control conditions for a valid quantitative run have not been met.

In case of an invalid quantitative run, repeat the PCR or extract DNA from the original samples again if no DNA is left over.

Qualitative analysis

A summary of results interpretation is shown in Table 6.

Table 6. Summary of results interpretation

Sample ID	Detection channel		Result interpretation
	Cycling Green	Cycling Yellow	
A	POSITIVE	POSITIVE*	JCV-specific DNA detected.
B	NEGATIVE	POSITIVE	JCV-specific DNA not detected. Sample does not contain detectable amounts of JCV-specific DNA.
C	NEGATIVE	NEGATIVE	PCR inhibition or reagent failure. Repeat procedure using original sample or collect and test a new sample.

* Detection of the Internal Control in the Cycling Yellow detection channel is not required for positive results in the Cycling Green detection channel. A high JCV load in the sample can lead to reduced or absent Internal Control signals.

Quantitative analysis

The *artus* JCV RG PCR Kit contains 4 Quantification Standards (QS). To generate a standard curve for quantitative analysis, these have to be defined as standards with appropriate concentrations (see Table 1, page 11). A standard curve for quantitative analysis can be generated using standards of known concentrations.

$$C_T = m \log(N_0) + b$$

- C_T = Threshold cycle
- m = Slope
- N_0 = Initial concentration
- b = Intercept

The concentrations of positive samples of unknown concentration can be derived from the standard curve (Figure 9).

$$N_0 = 10^{(C_T - b)/m}$$

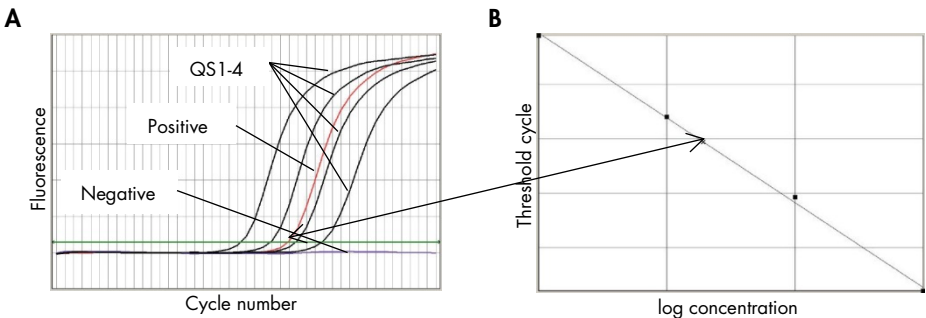


Figure 9. Quantification Standards, a positive and a negative sample displayed in (A) an amplification plot and (B) standard curve analysis.

Note: The concentration of the sample is displayed in IU/ μ l and refers to the concentration of viral DNA in the eluate.

Use the following formula to determine the viral load of the original sample.

$$\text{Viral load (sample) [IU/ml]} = \frac{\text{Volume (eluate) [\mu l]} \times \text{viral load (eluate) [IU/\mu l]}}{\text{Sample input [ml]}}$$

Limitations

- Use of this product is limited to personnel specially instructed and trained in the techniques of real-time PCR and in vitro diagnostic procedures.
- Good laboratory practice is essential for proper performance of this assay.
- Take extreme care to preserve the purity of the components of the kit and reaction setups. Closely monitor all reagents for impurities and contamination. Discard any reagents suspected of contamination.
- Appropriate specimen collection, transport, storage, and processing procedures are required for optimal performance of this assay.
- Do not use this assay directly on the specimen. Perform the applicable nucleic acid extraction procedures prior to using this assay.
- The presence of PCR inhibitors may cause false-negative or invalid results.
- Potential mutations within the target regions of the JCV genome covered by the primers and/or probes used in the kit may result in failure to detect the presence of the pathogen.
- As with any diagnostic test, interpret the results obtained using the *artus* JCV RG PCR Kit in consideration of all clinical and laboratory findings.

Quality Control

Each lot of *artus* JCV RG PCR Kit is tested against predetermined specifications to ensure consistent product quality.

Performance Characteristics

Performance evaluation of the *artus* JCV RG PCR Kit was carried out using quantified JCV-specific DNA isolated from JCV strain MAD-1 (ATCC® number: 45027).

Analytical sensitivity

The analytical sensitivity of the *artus* JCV RG PCR Kit is defined as the concentration (copies per μ l) of the eluate) of JCV-specific DNA that can be detected with a positivity rate of $\geq 95\%$. The analytical sensitivity was determined by analysis of a dilution series of quantified JCV DNA (Table 7).

Table 7. PCR results used to calculate the analytical sensitivity of the *artus* JCV RG PCR Kit

Input concentration (copies/ μ l)	Number of replicates	Number of positives	Hit rate (%)
3.2	18	18	100
1.0	18	18	100
0.3	18	18	100
0.1	18	18	100
0.03	18	16	89
0.01	18	13	72
0.003	18	15	83
0.001	36	29	81
0.0003	18	6	33
0.0001	18	3	17

The analytical sensitivity of the *artus* JCV RG PCR Kit, determined by probit analysis, for detection of JCV-specific DNA is 0.033 copies/ μ l (95% confidence interval [CI]: 0.013–0.175 copies/ μ l).

Analytical specificity

The analytical specificity of the *artus* JCV RG PCR Kit is ensured by careful selection of the oligonucleotides (primers and probes). The oligonucleotides are checked by sequence comparison analysis against publicly available sequences to ensure that all relevant JCV genotypes are detected.

The analytical specificity of the *artus* JCV RG PCR Kit was evaluated by testing a panel of genomic DNA/RNA extracted from other pathogens (Table 8).

Table 8. Organisms tested to demonstrate the analytical specificity of the *artus* JCV RG PCR Kit

Organism	Detection channel	
	Cycling Green (JCV)	Cycling Yellow (IC)
BK virus	Negative	Valid
Simian virus 40	Negative	Valid
Herpes simplex virus 1	Negative	Valid
Herpes simplex virus 2	Negative	Valid
Varicella-zoster virus	Negative	Valid
Epstein-Barr Virus	Negative	Valid
Cytomegalovirus	Negative	Valid
Human herpesvirus 6 (A, B)	Negative	Valid
Human herpesvirus 8	Negative	Valid
Hepatitis A virus	Negative	Valid
Hepatitis B virus	Negative	Valid
Hepatitis C virus	Negative	Valid
Human immunodeficiency virus 1	Negative	Valid
Parvovirus B19	Negative	Valid

The *artus* JCV RG PCR Kit did not cross-react with any of the specified organisms.

Linear range

The linear range of the *artus* JCV RG PCR Kit was evaluated by analyzing a logarithmic dilution series of DNA from JC virus strain MAD-1 using concentrations ranging from 1×10^9 to 10^{-1} copies/ μ l. At least 6 replicates per dilution were analyzed.

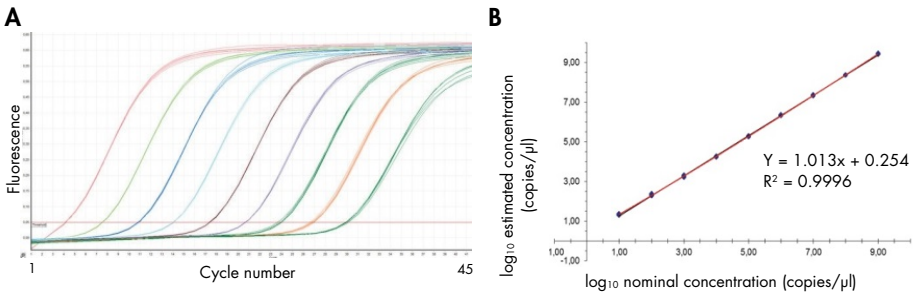


Figure 10. Amplification curve (A) and linear regression analysis (B) of a dilution series of DNA from JC virus strain MAD-1.

The linear range of the *artus* JCV RG PCR Kit extends over an interval of at least 8 orders of magnitude for JCV-specific DNA.

Precision

The precision of the *artus* JCV RG PCR Kit was determined as intra-assay variability (variability within one experiment), inter-assay variability (variability between different experiments), and inter-lot variability (variability between different production lots).

Variability data are expressed in terms of standard deviation, variance and coefficient of variation. The data are based on quantification analysis of defined concentrations of JCV-specific DNA and on threshold cycle (C_T) values in terms of the Internal Control (Tables 9 and 10). At least 6 replicates per sample were analyzed for intra-assay, inter-assay and inter-lot variability. Total variance was calculated by combining the 3 analyses.

Table 9. Precision data for the *artus* JCV RG PCR Kit (high positive sample)

JCV-specific system	Average conc. (copies/μl)	Standard deviation	Variance	Coefficient of variation (%)
Intra-assay variability	66.98	8.61	74.19	12.90
Inter-assay variability	72.93	8.85	78.30	12.13
Inter-lot variability	64.48	6.66	44.39	10.33
Total variance	69.27	9.02	81.39	13.02

Table 10. Precision data for the *artus* JCV RG PCR Kit (low positive sample)

JCV-specific system	Average conc. (copies/μl)	Standard deviation	Variance	Coefficient of variation (%)
Intra-assay variability	1.00	0.50	0.25	49.40
Inter-assay variability	1.08	0.43	0.19	40.00
Inter-lot variability	1.21	0.43	0.18	35.10
Total variance	1.14	0.45	0.20	39.10

Table 11. Precision data for the Internal Control of the *artus* JCV RG PCR Kit

Internal Control	Average threshold cycle (C _t)	Standard deviation	Variance	Coefficient of variation (%)
Intra-assay variability	22.02	0.16	0.03	0.74
Inter-assay variability	22.17	0.21	0.05	0.96
Inter-lot variability	22.26	0.28	0.08	1.25
Total variance	22.28	0.24	0.06	1.08

Repeatability







Specificity, sensitivity, and accuracy of quantification of the *artus* JCV RG PCR Kit were evaluated by analyzing established proficiency panels for JCV. To ensure repeatability of the *artus* JCV RG PCR Kit, specificity and sensitivity are evaluated by analyzing established proficiency panels for JCV as well as characterized diagnostic samples on a regular basis (Table 12).





Table 12. Results of the analysis of a proficiency panel for JC virus (QCMD) using the *artus* JC virus RG PCR Kit

Proficiency panel			<i>artus</i> JC virus RG PCR Kit		
Sample ID	Sample content	Expected conc. (copies/ml)	Result	Detected conc. of JC virus (copies/ml)	Internal Control
JCDNA14-01	JC virus type 1A	287	Positive	318	Valid
JCDNA14-02	JC virus type 1A	18,281	Positive	34,800	Valid
JCDNA14-03	JC virus type 3A	140	Positive	72	Valid
JCDNA14-04	JC virus 1A and BK virus 1b-2	977	Positive	948	Valid
JCDNA14-05	JC virus type 2B	1119	Positive	237	Valid
JCDNA14-06	BK virus type 1b-2	–	Negative	–	Valid
JCDNA14-07	JC virus type 3A	497	Positive	219	Valid
JCDNA14-08	JC virus type 1A	2178	Positive	4020	Valid
JCDNA14-09	JC virus negative	–	Negative	–	Valid
JCDNA14-10	JC virus type 2B	199	Positive	122	Valid

Symbols

The symbols in the following table are used in these instructions for use.

Symbol	Symbol definition
 96	Contains sufficient for 96 tests
	In vitro diagnostic medical device
	Catalog number
	Lot number
	Temperature limitation
	Manufacturer

Symbol	Symbol definition
	Use by
	Material number
	Global Trade Item Number
	Consult instructions for use

Troubleshooting Guide

The scientists in QIAGEN Technical Services are always happy to answer any questions you may have about either the information and/or protocols in this handbook or sample and assay technologies (for contact information, visit www.qiagen.com).

Ordering Information

Product	Contents	Cat. no.
<i>artus</i> JCV RG PCR Kit (96)	For 96 reactions: Master A, Master B, 4 Quantification Standards, Internal Control, H ₂ O (PCR grade water)	4532265
QIAamp DNA Mini Kit (50)	For 50 DNA preps: 50 QIAamp Mini Spin Columns, Proteinase K, Reagents, Buffers, Collection Tubes (2 ml)	51304
QIAamp DNA Mini Kit (250)	For 250 DNA preps: 250 QIAamp Mini Spin Columns, Proteinase K, Reagents, Buffers, Collection Tubes (2 ml)	51306
Rotor-Gene Q and accessories		
Rotor-Gene Q MDx 5plex System	Real-time PCR cyclers with 5 channels (green, yellow, orange, red, crimson), laptop computer, software, accessories: includes 1-year warranty on parts and labor, installation and training	9002023
Rotor-Gene Q MDx 5plex Platform	Real-time PCR cyclers with 5 channels (green, yellow, orange, red, crimson), laptop computer, software, accessories: includes 1-year warranty on parts and labor, installation and training not included	9002022

Product	Contents	Cat. no.
Rotor-Gene Q 5plex Priority Package Plus	Real-time PCR cyclers with 5 channels (green, yellow, orange, red, crimson), laptop computer, software, accessories: includes Priority Package with software, installation, training, 3-year warranty on parts and labor, and 3 preventive maintenance visits	9001866
Rotor-Gene Q 5plex Priority Package	Real-time PCR cyclers with 5 channels (green, yellow, orange, red, crimson), laptop computer, software, accessories: includes Priority Package with software, installation, training, 2-year warranty on parts and labor, and 2 preventive maintenance visits	9001865
Rotor-Gene Q 5plex System	Real-time PCR cyclers with 5 channels (green, yellow, orange, red, crimson), laptop computer, software, accessories: includes 1-year warranty on parts and labor, installation and training	9001640
Rotor-Gene Q 5plex Platform	Real-time PCR cyclers with 5 channels (green, yellow, orange, red, crimson), laptop computer, software, accessories: includes 1-year warranty on parts and labor, installation and training not included	9001570

Product	Contents	Cat. no.
Rotor-Gene Q 6plex Priority Package Plus	Real-time PCR instrument with 6 channels (blue, green, yellow, orange, red, crimson), including laptop computer, software, accessories: includes Priority Package with software, installation, training, 3-year warranty on parts and labor, and 3 preventive maintenance visits	9001870
Rotor-Gene Q 6plex Priority Package	Real-time PCR instrument with 6 channels (blue, green, yellow, orange, red, crimson), including laptop computer, software, accessories: includes Priority Package with software, installation, training, 2-year warranty on parts and labor, and 2 preventive maintenance visits	9001869
Rotor-Gene Q 6plex System	Real-time PCR instrument with 6 channels (blue, green, yellow, orange, red, crimson), including laptop computer, software, accessories: includes 1-year warranty on parts and labor, installation and training	9001660
Rotor-Gene Q 6plex Platform	Real-time PCR instrument with 6 channels (blue, green, yellow, orange, red, crimson), including laptop computer, software, accessories: includes 1-year warranty on parts and labor, installation and training not included	9001590

Product	Contents	Cat. no.
Loading Block 72 x 0.1 ml Tubes	Aluminum block for manual reaction setup with a single-channel pipet in 72 x 0.1 ml tubes	9018901
Strip Tubes and Caps, 0.1 ml (250)	250 strips of 4 tubes and caps for 1000 reactions of 10–50 µl	981103
Strip Tubes and Caps, 0.1 ml (2500)	10 x 250 strips of 4 tubes and caps for 10,000 reactions of 10–50 µl	981106

Notes

Notes

Handbook Revision History

Date	Changes
R2	Changed assay reporting units from copies to International Units (IU).
R3	Minor error corrections: update to release date.

Limited License Agreement for *artus JCV RG PCR Kit*

Use of this product signifies the agreement of any purchaser or user of the product to the following terms:

1. The product may be used solely in accordance with the protocols provided with the product and this handbook and for use with components contained in the kit only. QIAGEN grants no license under any of its intellectual property to use or incorporate the enclosed components of this kit with any components not included within this kit except as described in the protocols provided with the product, this handbook, and additional protocols available at www.qiagen.com. Some of these additional protocols have been provided by QIAGEN users for QIAGEN users. These protocols have not been thoroughly tested or optimized by QIAGEN. QIAGEN neither guarantees them nor warrants that they do not infringe the rights of third-parties.
2. Other than expressly stated licenses, QIAGEN makes no warranty that this kit and/or its use(s) do not infringe the rights of third-parties.
3. This kit and its components are licensed for one-time use and may not be reused, refurbished, or resold.
4. QIAGEN specifically disclaims any other licenses, expressed or implied other than those expressly stated.
5. The purchaser and user of the kit agree not to take or permit anyone else to take any steps that could lead to or facilitate any acts prohibited above. QIAGEN may enforce the prohibitions of this Limited License Agreement in any Court, and shall recover all its investigative and Court costs, including attorney fees, in any action to enforce this Limited License Agreement or any of its intellectual property rights relating to the kit and/or its components.

For updated license terms, see www.qiagen.com.

The purchase of this product allows the purchaser to use it for the performance of diagnostic services for human in vitro diagnostics. No general patent or other license of any kind other than this specific right of use from purchase is granted hereby.

Trademarks: QIAGEN®, Sample to Insight®, QIAamp®, *artus*®, Rotor-Gene® (QIAGEN Group); ATCC® (American Type Culture Collection Corporation); FAM™, JOE™ (Thermo Fisher Scientific or its subsidiaries).

1121142 03/2020 HB-2017-004 © 2020 altona Diagnostics GmbH, all rights reserved.

Ordering www.qiagen.com/contact | Technical Support support.qiagen.com | Website www.qiagen.com