

REF **900202 NeuMoDx™ HCV External Controls**

R only

CAUTION: For US Export Only

IVD For *in vitro* diagnostic use with the NeuMoDx 288 and NeuMoDx 96 Molecular System



For insert updates, go to: www.qiaagen.com/neumodx-ifu

For detailed instructions, refer to the NeuMoDx 288 Molecular System Operator's Manual; p/n 40600108

For detailed instructions, refer to the NeuMoDx 96 Molecular System Operator's Manual; p/n 40600317

See also the NeuMoDx HCV Quant Test Strip Instructions for Use; p/n 40600140

INTENDED USE

The NeuMoDx HCV External Controls are a component of the NeuMoDx HCV Quant Assay, an *in vitro* diagnostic nucleic acid amplification test intended for the detection and quantitation of hepatitis C virus (HCV) RNA in human plasma. As implemented on the fully automated NeuMoDx 288 Molecular System or NeuMoDx 96 Molecular System (NeuMoDx System(s)), the NeuMoDx HCV External Controls are used to establish runtime validity required to execute the NeuMoDx HCV Quant Assay for accurate quantitation of HCV RNA in human plasma specimens.

SUMMARY AND EXPLANATION

The NeuMoDx HCV External Controls are provided in 15 paired sets of positive and negative control vials. One set of external controls is processed every 24 hours to establish runtime validity of the NeuMoDx HCV Quant Assay. The HCV target in the positive control is a non-infectious, replication-defective mammalian recombinant virus containing HCV genome sequences and diluted in Basematrix 53 Diluent (Basematrix) (Seracare Life Sciences, Milford, MA, USA). The negative HCV control consists of Basematrix only.

The NeuMoDx HCV Quant Assay combines automated RNA extraction, amplification, and detection by real-time reverse transcriptase PCR to enable the quantitative detection of HCV RNA in human plasma specimens. The NeuMoDx HCV Quant Assay includes an exogenous RNA Sample Process Control (SPC2) to help monitor for the presence of potential inhibitory substances and for NeuMoDx System or reagent failures that may be encountered during the extraction and amplification processes.

Clinical laboratories typically require that external controls be incorporated into routine testing protocols to assess test performance and ensure that the test procedures meet established quality control requirements. The NeuMoDx HCV External Controls are used to establish such routine run validity of the NeuMoDx HCV Quant Assay. Routine use of these controls enables the laboratories to monitor day-to-day variation and lot-to-lot performance of the NeuMoDx HCV Quant Assay reagents and can assist the lab in identifying errors prior to reporting of test results.

PRINCIPLES OF THE PROCEDURE

The NeuMoDx HCV External Controls are non-infectious samples formulated to mimic naturally occurring human plasma specimens. The encapsulated target material used in the positive control allows for the verification of efficacious nucleic acid extraction procedure. One set of controls is processed every 24 hours. Such routine processing of the NeuMoDx HCV External Controls enables the laboratories to ensure reliability of test results for human clinical specimens processed within the 24-hour validity period. The external controls are processed in a manner identical to the processing of the human clinical specimens intended for quantitative HCV testing.

Expected results for both these external controls are incorporated into the Control Validity algorithm included in the NeuMoDx System software. Upon successful processing of the external controls, the system software automatically records the validity for a period of 24 hours. The system software automatically alerts the user to process the external controls when control validity period has expired.



REAGENTS / CONSUMABLES

Material Provided

REF	Contents	Tests per unit	Total tests per kit
900202	NeuMoDx HCV External Controls <i>Single use sets of HCV Positive and Negative Controls to establish daily validity of NeuMoDx HCV Quant Assay (1 vial of each control = 1 set)</i>	1 set	15

Materials Required but Not Provided (Available Separately from NeuMoDx)

REF	Contents
300300	NeuMoDx HCV Quant Test Strip <i>Dried PCR reagents containing HCV and SPC2 specific TaqMan® probes and primers</i>
100200	NeuMoDx Extraction Plate <i>Dried paramagnetic particles, lytic enzyme, and sample process controls</i>
800202	NeuMoDx HCV Calibrators <i>Single use sets of HCV High and Low Calibrators to establish validity of standard curve</i>
400600	NeuMoDx Lysis Buffer 3
400100	NeuMoDx Wash Reagent
400200	NeuMoDx Release Reagent
100100	NeuMoDx Cartridge
235903	Hamilton CO-RE / CO-RE II Tips (300 µL) with Filters
235905	Hamilton CO-RE / CO-RE II Tips (1000 µL) with Filters

Instrumentation Required

NeuMoDx 288 Molecular System [REF 500100] or NeuMoDx 96 Molecular System [REF 500200]



WARNINGS AND PRECAUTIONS

- The NeuMoDx HCV External Controls are for *in vitro* diagnostic use only with the NeuMoDx HCV Quant Test Strip as implemented on the NeuMoDx System.
- Do not use the NeuMoDx HCV External Controls after the listed expiration date.
- Do not use the NeuMoDx HCV External Controls if the packaging is damaged or the contents are not frozen upon arrival.
- Always handle specimens as if they are infectious and in accordance with safe laboratory procedures such as those described in *Biosafety in Microbiological and Biomedical Laboratories*¹ and in CLSI Document M29-A4.²
- Do not pipette by mouth. Do not smoke, drink, or eat in areas where specimens or reagents are being handled.
- Dispose of unused reagents and waste in accordance with country, federal, provincial, state and local regulations.
- Clean, powder-free, nitrile gloves should be worn when handling all NeuMoDx reagents and consumables.
- Wash hands thoroughly after performing the test.
- Safety Data Sheets (SDS) are provided for each reagent (as applicable) at www.giagen.com/neumodx-ifu
- Do not reuse.



PRODUCT STORAGE, HANDLING AND STABILITY

- The NeuMoDx HCV External Controls are shipped with dry ice to maintain a frozen state; do not use if contents are not frozen upon receipt.
- It is recommended that the NeuMoDx HCV External Controls be stored at -15 °C to -20 °C to ensure stability. NeuMoDx HCV External Controls' shelf life is 24 months from the date of manufacture.
- Control vials are intended for single use only. Thawed external controls may be stored at 4 °C for no longer than 24 hours.
- Refreezing after a first thaw is not recommended.
- Although the NeuMoDx HCV External Controls are non-infectious, any unused material should be discarded after use as biohazard waste to reduce risk of contamination by the target nucleic acid contained.
- Discard any controls that appear cloudy or contain large precipitates after thawing.

INSTRUCTIONS FOR USE

1. One set of external controls needs to be processed every 24 hours throughout testing with the NeuMoDx HCV Quant Assay. If a set of valid test controls does not exist, the NeuMoDx System software will prompt the user for these controls to be processed before sample results can be reported.
2. If external controls are required, process the controls (1 positive control and 1 negative control):

NeuMoDx HCV External Control	Label Color Scheme
Positive Control (HCVPC)	Red
Negative Control (HCVNC)	Black

3. Retrieve the set of NeuMoDx HCV External Controls from freezer and allow the vials to set at room temperature (15-30 °C) until completely thawed.
4. Vortex gently to ensure homogeneity.
5. Load the control vials into a standard 32-tube Specimen Tube Carrier, and ensure caps are removed from all tubes.
6. Place the Specimen Tube Carrier on the Autoloader shelf and use the touchscreen to load carrier into the NeuMoDx System.
7. The NeuMoDx System will recognize the barcode and start processing the specimen tubes unless reagents or consumables required for testing are not available.
8. Validity of these external controls will be assessed by the NeuMoDx System based on the expected results.

NeuMoDx HCV External Control	HCV Result	SPC2 Result
Positive Control (HCVPC)	HCV POSITIVE	N/A
Negative Control (HCVNC)	HCV NEGATIVE	SPC2 Positive

9. Discrepant result handling for external controls should be performed as follows:
 - a) A Positive test result reported for a negative control sample indicates a specimen contamination problem.
 - b) A negative result reported for a positive control sample may indicate there is a reagent or instrument related problem.
 - c) In either of the above instances, or in the event of an indeterminate (IND) result, repeat the failed control with freshly thawed vial(s) of the control(s) failing the validity test.
 - d) If the Positive external control continues to report a Negative result, contact NeuMoDx technical support.
 - e) If the Negative external control continues to report a Positive result, attempt to eliminate all sources of potential contamination, including replacing all reagents and repeat the run before contacting NeuMoDx customer service.

LIMITATIONS

1. The NeuMoDx HCV External Controls can only be used in conjunction with NeuMoDx HCV Quant Test Strip on the NeuMoDx Systems.
2. A valid calibration of the NeuMoDx HCV Quant Test Strip using NeuMoDx HCV External Calibrators is required *before* the NeuMoDx HCV External Controls can be processed.
3. Erroneous results could occur from improper handling, storage, or other technical error.
4. Operation of the NeuMoDx System is limited to use by personnel trained on the use of the NeuMoDx System.

REFERENCES

1. Biosafety in Microbiological and Biomedical Laboratories, 5th edition. HHS Publication No. (CDC) 21-1112, Revised December 2009.
2. Clinical and Laboratory Standards Institute (CLSI). Protection of Laboratory Workers from Occupationally Acquired Infections; Approved Guideline – Fourth Edition. CLSI document M29-A4; May 2014.





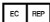


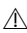
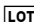



TRADEMARKS

NeuMoDx is a trademark of NeuMoDx Molecular, Inc.

TaqMan® is a registered trademark of Roche Molecular Systems, Inc.

All other product names, trademarks, and registered trademarks that may appear in this document are property of their respective owners.

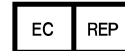
SYMBOL KEY

R only	Prescription use only		Temperature limit
	Manufacturer		Do not re-use
	<i>In vitro</i> diagnostic medical device		Contains sufficient for <n> tests
	Authorized representative in the European Community		Consult instructions for use
	Catalog number		Caution
	Batch code		Biological risks
	Use-by date		CE Mark

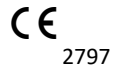


NeuMoDx Molecular, Inc.
1250 Eisenhower Place
Ann Arbor, MI 48108, USA

Sponsor (AUS):
QIAGEN Pty Ltd
Level 2 Chadstone Place
1341 Dandenong Rd
Chadstone VIC 3148
Australia



Emergo Europe B.V.
Westervoortsedijk 60
6827 AT Arnhem
The Netherlands



Technical support/Vigilance reporting: support.qiagen.com

Patent: www.neumodx.com/patents