

August 2015

---

# GeneRead™ Library Quant Array Handbook

For reliable quantification of Ion Torrent™ or  
Illumina® libraries using real-time PCR



---

Sample & Assay Technologies

## QIAGEN Sample and Assay Technologies

QIAGEN is the leading provider of innovative sample and assay technologies, enabling the isolation and detection of contents of any biological sample. Our advanced, high-quality products and services ensure success from sample to result.

### QIAGEN sets standards in:

- Purification of DNA, RNA and proteins
- Nucleic acid and protein assays
- microRNA research and RNAi
- Automation of sample and assay technologies

Our mission is to enable you to achieve outstanding success and breakthroughs. For more information, visit [www.qiagen.com](http://www.qiagen.com).

# Contents

<b>Kit Contents</b>	<b>4</b>
<b>Storage</b>	<b>5</b>
<b>Quality Control</b>	<b>5</b>
<b>Intended Use</b>	<b>5</b>
<b>Technical Assistance</b>	<b>5</b>
<b>Introduction</b>	<b>7</b>
Principle and Workflow	7
<b>Equipment and Reagents to Be Supplied by User</b>	<b>10</b>
<b>Important Notes</b>	<b>11</b>
DNA contamination	11
General considerations	11
<b>Protocol 1: Real-Time PCR for GeneRead Library Quant Array for Ion Torrent or Illumina (Formats A, C, D, E, F, G)</b>	<b>13</b>
<b>Protocol 2: Real-Time PCR for GeneRead Library Quant Array for Ion Torrent or Illumina (Format R)</b>	<b>18</b>
<b>Data Analysis</b>	<b>22</b>
Obtaining raw threshold cycle ( $C_T$ ) values	22
Exporting $C_T$ values	22
Excel-based data analysis file	22
<b>Troubleshooting Guide</b>	<b>24</b>
<b>References</b>	<b>25</b>
<b>Ordering Information</b>	<b>26</b>

## Kit Contents

<b>GeneRead Library Quant Array</b>							
<b>Catalog no.</b>	<b>180611 (Varies depending on format)</b>						
<b>Format</b>	<b>A</b>	<b>C</b>	<b>D</b>	<b>E</b>	<b>F</b>	<b>G</b>	<b>R</b>
Plate containing dried assays	96-well (2)	96-well (2)	96-well (2)	384-well (1)	96-well (2)	384-well (1)	100-well (2)
Optical Thin-Wall 8-Cap Strips (12 per plate)	12 x 2	–	12 x 2	–	–	–	–
Optical adhesive film	–	1 x 2	–	1	1 x 2	1	–
Rotor-Disc Heat-Sealing Film	–	–	–	–	–	–	1 x 2
Handbook				1			

### Array formats for use with real-time cyclers

<b>Format</b>	<b>Suitable real-time cyclers</b>	<b>Plate</b>
A	Applied Biosystems® 5700, 7000, 7300, 7500 Standard, 7700, 7900HT Standard; Bio-Rad® iCycler®, iQ™ 5, MyiQ™, MyiQ2™, Bio-Rad/MJ Research Chromo4™; Eppendorf® Mastercycler® ep realplex 2, 2s, 4, 4s; Stratagene® Mx3005P®, Mx3000P®	96-well
C	Applied Biosystems 7500 FAST, 7900HT FAST, StepOnePlus™	96-well
D	Bio-Rad CFX96™, Bio-Rad/MJ Research Opticon 2®; Stratagene Mx4000®	96-well
E	Applied Biosystems 7900HT (384-well block); Bio-Rad CFX384™	384-well
F	Roche® LightCycler® 480 II (96-well block)	96-well
G	Roche LightCycler 480 II (384-well block)	384-well
R	QIAGEN Rotor-Gene® Q and Rotor-Gene 6000	100-well

## Storage

The GeneRead Library Quant Array is shipped at ambient temperature, on ice, or on dry ice depending on the destination and accompanying products. Upon receipt, store at  $-20^{\circ}\text{C}$ . If stored under these conditions, GeneRead Library Quant Array is stable for 6 months after receipt.

## Quality Control

In accordance with QIAGEN's ISO-certified Quality Management System, each lot of GeneRead Library Quant Array is tested against predetermined specifications to ensure consistent product quality.

## Intended Use

The GeneRead Library Quant Array is intended for molecular biology applications. This product is not intended for the diagnosis, prevention or treatment of a disease.

All due care and attention should be exercised in the handling of the products. We recommend all users of QIAGEN products to adhere to the NIH guidelines that have been developed for recombinant DNA experiments, or to other applicable guidelines.

## Technical Assistance

At QIAGEN, we pride ourselves on the quality and availability of our technical support. Our Technical Service Departments are staffed by experienced scientists with extensive practical and theoretical expertise in sample and assay technologies and the use of QIAGEN products. If you have any questions or experience any difficulties regarding GeneRead Library Quant Array, or QIAGEN products in general, please do not hesitate to contact us.

QIAGEN customers are a major source of information regarding advanced or specialized uses of our products. This information is helpful to other scientists as well as to the researchers at QIAGEN. We therefore encourage you to contact us if you have any suggestions about product performance or new applications and techniques.

For technical assistance and more information, please see our Technical Support Center at [www.qiagen.com/Support](http://www.qiagen.com/Support) or call one of the QIAGEN Technical Service Departments or local distributors (see back cover or visit [www.qiagen.com](http://www.qiagen.com)).

## Safety Information

When working with chemicals, always wear a suitable lab coat, disposable gloves and protective goggles. For more information, please consult the appropriate safety data sheets (SDSs). These are available online in convenient and compact PDF format at [www.qiagen.com/safety](http://www.qiagen.com/safety) where you can find, view and print the SDS for each QIAGEN kit and kit component.

24-hour emergency information

Chemical emergency or accident assistance is available 24 hours a day from:

CHEMTREC

USA & Canada ■ Tel: 1-800-424-9300

Outside USA & Canada ■ Tel: +1-703-527-3887 (collect calls accepted)

## Introduction

One of the most important factors in a next-generation sequencing experiment is the accurate quantification of the prepared library. The accurate quantification of amplifiable library molecules is essential for ensuring optimal quality reads and efficient data generation. Underestimation of amplifiable library molecules leads to mixed signals and non-resolvable data; conversely, overestimation results in poor yield of template-carrying beads (Ion Torrent platform) or clusters (Illumina platform) and reduced usage of sequencing capacity.

GeneRead Library Quant Array uses real-time PCR to quantify NGS libraries. GeneRead Library Quant Array specifically quantifies DNA molecules with adaptors at both ends, which are the only amplifiable molecules during emulsion PCR (Ion Torrent platform) or bridge PCR (Illumina platform) and therefore provides highly accurate quantification of amplifiable library molecules. The high sensitivity of real-time PCR allows quantification of libraries with very low concentrations, even below the detection threshold of conventional spectrophotometric methods.

The GeneRead Library Quant Array is optimized with GeneRead qPCR Mastermixes to provide superior sensitivity and wide linear dynamic ranges. It can be easily automated for high-throughput applications.

## Principle and Workflow

The GeneRead Library Quant Array for Ion Torrent contains five predispensed, sequential 10-fold dilutions of Ion Torrent DNA Standard mixed with a PCR primer assay in triplicates, as well as PCR primer assays in the remaining wells of a 96-well, 384-well or 100-well PCR plate. The predispensed serially diluted DNA standards and PCR primer assay provide a highly convenient method for quantification of library input.

The Ion Torrent DNA Standard harbors a 130 bp target flanked with the "A" and "trP1" adaptor sequences used in Ion Torrent sequencing libraries. The Ion Torrent DNA Standard generates a 183 bp amplicon using the provided PCR assay. This PCR assay is experimentally verified to specifically amplify Ion Torrent library molecules flanked with the "A" and "trP1" adaptor sequences. The primer sequences used in the PCR assay are as follows:

Ion Torrent forward: 5'-CCA TCT CAT CCC TGC GTG TC-3'

Ion Torrent reverse: 5'-CCT CTC TAT GGG CAG TCG GTG AT-3'

The GeneRead Library Quant Array for Illumina contains five predispensed, sequential 10-fold dilutions of Illumina DNA Standard mixed with a PCR primer assay in triplicates, as well as PCR primer assays in the remaining wells of a 96-well, 384-well or 100-well PCR plate. The predispensed, serially diluted DNA standards and PCR primer assay provide a highly convenient method for quantification of library input.


The Illumina DNA Standard harbors a target flanked by the "P5" and "P7" primer sequences used in Illumina sequencing libraries. The Illumina DNA standard generates a 426 bp amplicon using the provided PCR assay. This PCR assay is experimentally verified to specifically amplify Illumina library molecules with the "P5" and "P7" sequences. The primer sequences used in the PCR assay are as follows:

Illumina forward: 5'-AAT GAT ACG GCG ACC ACC GA-3'

Illumina reverse: 5'-CAA GCA GAA GAC GGC ATA CGA-3'

The GeneRead Library Quant Array workflow (Figure 1) is simple. The procedure begins with two 10-fold dilutions of the sample library (the sample library must be diluted to a concentration within the range of the serially diluted standards). Next, appropriate mixtures of PCR reactions are prepared and aliquotted into the wells of a real-time PCR plate. PCR is performed and  $C_T$  values are exported to the provided Excel<sup>®</sup>-based data analysis file to calculate the Template Dilution Factor (for the Ion Torrent platform) or library concentration (for the Illumina platform).



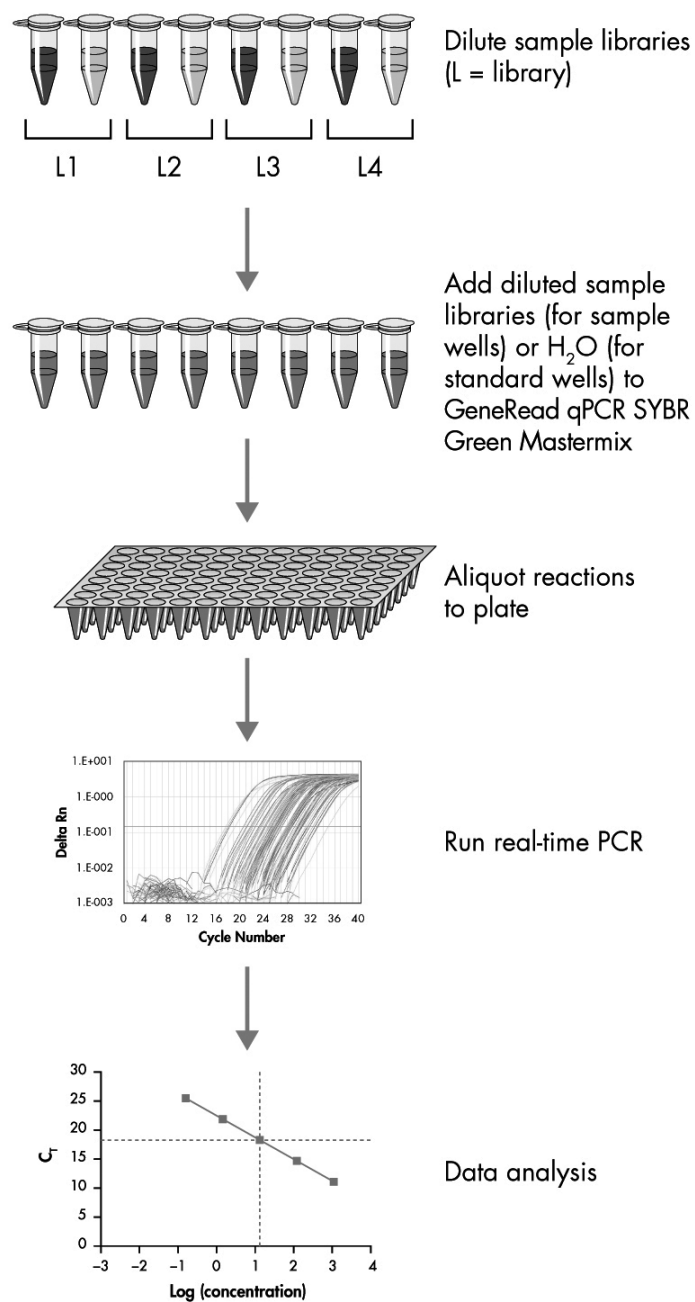

**Figure 1. GeneRead Library Quant Array workflow.**

## Equipment and Reagents to Be Supplied by User

When working with chemicals, always wear a suitable lab coat, disposable gloves and protective goggles. For more information, consult the appropriate safety data sheets (SDSs), available from the product supplier.

- Appropriate GeneRead qPCR SYBR® Green Mastermix (be sure to select the correct format for the PCR instrument):

Master mix	Cat. no.	Instrument
GeneRead qPCR SYBR Green ROX™ Mastermix	180840	All Applied Biosystems and Stratagene instruments, and Eppendorf Mastercycler ep <i>realplex</i> instruments with a ROX filter set
GeneRead qPCR SYBR Green Fluor Mastermix	180830	Bio-Rad iCycler, MyiQ, MyiQ2 and iQ 5
GeneRead qPCR SYBR Green Mastermix	180820	Bio-Rad models CFX96, CFX384, Bio-Rad/MJ Research Opticon 2, and Bio-Rad/MJ Research Chromo4, Roche LightCycler 480 II (96- and 384-well)
GeneRead qPCR SYBR Green ROX FAST Mastermix	180850	Rotor-Gene Q and Rotor-Gene 6000

- Real-time PCR instrument
- Calibrated single- and multi-channel pipets
- RNase-/DNase-free pipet tips and tubes
- RNase-/DNase-free 200 µl regular PCR tubes, or 8-tube or 12-tube strips
- Molecular biology grade RNase- and DNase-free water
- Dilution buffer (10 mM Tris-Cl, pH 7.8 or QIAGEN EB buffer).
- Rotor-Disc® heat sealer (if using RGQ).

# Important Notes

## DNA contamination

For reliable results, it is very important to prevent contamination with foreign DNA. Even very small amounts of foreign DNA can artificially inflate SYBR Green signals, yielding false positive results. The most common source of contamination in PCR reagents comes from the products of previous PCR experiments in the working area. To minimize contamination, follow the recommendations below:

- Wear gloves throughout the entire procedure.
- Use only fresh PCR-grade reagents and labware.
- Physically separate the workspace for PCR setup and post-PCR work.
- Before setting up an experiment, decontaminate the PCR workspace and labware (pipet barrels, tube racks, etc.) with 10% bleach and UV light. Preferentially set up reactions in a PCR workstation.
- Close all tubes containing PCR products as soon as possible after use.
- Treat any labware (tips or tubes) containing PCR products or other DNA with 10% bleach before discarding.

## General considerations:

- For accuracy and precision, ensure that micropipettors are calibrated before beginning the protocol. Be sure not to introduce bubbles into the wells when pipetting.
- Do not use DEPC-treated water. Use high-quality, nuclease-free water.
- If precipitates are present in the master mix tubes, warm the reagents at 42°C for 1 min and vortex briefly to dissolve. Repeat if necessary.
- The GeneRead Library Quant Arrays have been tested on the following cyclers:

<b>Cycler</b>	<b>Plate</b>
Vii7A/Stratagene	A
ABI7500	C
ABI7900	E
Roche-LightCycler II	F
RGQ	Rotor-Disc 100

- Refer to tables 3, 4, 5 and 8 in this handbook for cycler-specific cycling programs.
- Instrument setup protocols can be found at:  
<https://www.qiagen.com/us/shop/assay-technologies/real-time-pcr-and-rt-pcr-reagents/rt2-profiler-pcr-arrays/#resources>

# Protocol 1: Real-Time PCR for GeneRead Library Quant Array for Ion Torrent or Illumina (Formats A, C, D, E, F, G)

## Important points before starting

- Ensure that the GeneRead qPCR SYBR Green Mastermix (see page 10) and the plate format are suitable for your real-time cycler (see page 4). An incorrect format will not fit the real-time cycler properly and may damage the real-time cycler.
- If you are dealing with a new library prep without prior measurement of the final concentration, we recommend checking the concentration using the Agilent® BioAnalyzer or QIAxcel® Advanced and diluting the library to about 1 nM before proceeding with the protocol given below.

## Procedure

### Preparing sample library dilution

1. Prepare a starting 1:20 dilution of the sample library in dilution buffer (2 µl sample library + 38 µl dilution buffer).
2. Prepare two working dilutions of sample library, Dilution 1 and Dilution 2, as indicated in Table 1.

**Table 1. Sample library dilution**

Dilution	Library	Dilution buffer
Dilution 1 (1:2000)	2 µl of 1:20	198 µl
Dilution 2 (1:20000)	5 µl of 1:2000	45 µl

**Note:** Always use freshly made dilutions.

### Setting up the PCR

3. Briefly centrifuge the GeneRead qPCR SYBR Green Mastermix (10–15 s) to bring the contents to the bottom of the tube.

**Note:** As the GeneRead qPCR SYBR Green Mastermix contains HotStarTaq® DNA Polymerase, which is active only after heat activation, reactions can be prepared at room temperature (15–25°C).

4. Make PCR mix for each sample dilution (for triplicates on A, C, D, F format plate and quadruplicates on E, G format plate) or standards / NTC wells according to Table 2.

**Table 2. PCR mix for standard or sample libraries\***

Component	For 96-well plate (A,C,D,F)		For 384-well plate (E,G)	
	For standard and NTC  (18 wells)	For sample library  Dilution 1 and 2 (3 wells for each dilution)	For standard and NTC	For sample library  Dilution 1 and 2 (4 wells for each dilution)
RNase-/DNase-free water	240 µl	32 µl	150 µl	19 µl
GeneRead qPCR SYBR Green Mastermix	240 µl	42 µl	150 µl	25 µl
Template (sample library)	–	10 µl	–	6 µl
<b>Final volume</b>	<b>480 µl</b>	<b>84 µl</b>	<b>300 µl</b>	<b>50 µl</b>

\* For multiple wells, an excess volume is provided to allow for pipetting errors. Perform pipetting steps as precisely as possible to ensure that each well receives the required volume.

5. Add 25 µl per well PCR mix to three wells for 96-well plate (see Figure 2 for layout) or 10 µl/per well PCR mix to four wells for 384-well plate (see Figure 3 for layout) for each sample dilution. Add 25 µl (96-well plate, Figure 2) or 10 µl (384-well plate, Figure 3) PCR reaction mix (with no template) into each well for standards and NTC.

The layouts of figures 2 and 3 must be followed as shown for appropriate data analysis.

Well	1	2	3	4	5	6	7	8	9	10	11	12
A	Standard 1	Standard 1	Standard 1	Library 2 dilution 1	Library 2 dilution 1	Library 2 dilution 1	Library 6 dilution 1	Library 6 dilution 1	Library 6 dilution 1	Library 10 dilution 1	Library 10 dilution 1	Library 10 dilution 1
B	Standard 2	Standard 2	Standard 2	Library 2 dilution 2	Library 2 dilution 2	Library 2 dilution 2	Library 6 dilution 2	Library 6 dilution 2	Library 6 dilution 2	Library 10 dilution 2	Library 10 dilution 2	Library 10 dilution 2
C	Standard 3	Standard 3	Standard 3	Library 3 dilution 1	Library 3 dilution 1	Library 3 dilution 1	Library 7 dilution 1	Library 7 dilution 1	Library 7 dilution 1	Library 11 dilution 1	Library 11 dilution 1	Library 11 dilution 1
D	Standard 4	Standard 4	Standard 4	Library 3 dilution 2	Library 3 dilution 2	Library 3 dilution 2	Library 7 dilution 2	Library 7 dilution 2	Library 7 dilution 2	Library 11 dilution 2	Library 11 dilution 2	Library 11 dilution 2
E	Standard 5	Standard 5	Standard 5	Library 4 dilution 1	Library 4 dilution 1	Library 4 dilution 1	Library 8 dilution 1	Library 8 dilution 1	Library 8 dilution 1	Library 12 dilution 1	Library 12 dilution 1	Library 12 dilution 1
F	NTC	NTC	NTC	Library 4 dilution 2	Library 4 dilution 2	Library 4 dilution 2	Library 8 dilution 2	Library 8 dilution 2	Library 8 dilution 2	Library 12 dilution 2	Library 12 dilution 2	Library 12 dilution 2
G	Library 1 dilution 1	Library 1 dilution 1	Library 1 dilution 1	Library 5 dilution 1	Library 5 dilution 1	Library 5 dilution 1	Library 9 dilution 1	Library 9 dilution 1	Library 9 dilution 1	Library 13 dilution 1	Library 13 dilution 1	Library 13 dilution 1
H	Library 1 dilution 2	Library 1 dilution 2	Library 1 dilution 2	Library 5 dilution 2	Library 5 dilution 2	Library 5 dilution 2	Library 9 dilution 2	Library 9 dilution 2	Library 9 dilution 2	Library 13 dilution 2	Library 13 dilution 2	Library 13 dilution 2

Figure 2. GeneRead Library Quant Array Layout (96-well plate).

Well	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24
A	standard1	Library 2 dilution 1	Library 6 dilution 1	Library 10 dilution 1	Library 14 dilution 1	Library 18 dilution 1	Library 22 dilution 1	Library 26 dilution 1	Library 30 dilution 1	Library 34 dilution 1	Library 38 dilution 1	Library 42 dilution 1												
B	standard2	Library 2 dilution 2	Library 6 dilution 2	Library 10 dilution 2	Library 14 dilution 2	Library 18 dilution 2	Library 22 dilution 2	Library 26 dilution 2	Library 30 dilution 2	Library 34 dilution 2	Library 38 dilution 2	Library 42 dilution 2												
C	standard3	Library 3 dilution 1	Library 7 dilution 1	Library 11 dilution 1	Library 15 dilution 1	Library 19 dilution 1	Library 23 dilution 1	Library 27 dilution 1	Library 31 dilution 1	Library 35 dilution 1	Library 39 dilution 1	Library 43 dilution 1												
D	standard4	Library 3 dilution 2	Library 7 dilution 2	Library 11 dilution 2	Library 15 dilution 2	Library 19 dilution 2	Library 23 dilution 2	Library 27 dilution 2	Library 31 dilution 2	Library 35 dilution 2	Library 39 dilution 2	Library 43 dilution 2												
E	standard5	Library 4 dilution 1	Library 8 dilution 1	Library 12 dilution 1	Library 16 dilution 1	Library 20 dilution 1	Library 24 dilution 1	Library 28 dilution 1	Library 32 dilution 1	Library 36 dilution 1	Library 40 dilution 1	Library 44 dilution 1												
F	NTC	Library 4 dilution 2	Library 8 dilution 2	Library 12 dilution 2	Library 16 dilution 2	Library 20 dilution 2	Library 24 dilution 2	Library 28 dilution 2	Library 32 dilution 2	Library 36 dilution 2	Library 40 dilution 2	Library 44 dilution 2												
G	Library 1 dilution 1	Library 5 dilution 1	Library 9 dilution 1	Library 13 dilution 1	Library 17 dilution 1	Library 21 dilution 1	Library 25 dilution 1	Library 29 dilution 1	Library 33 dilution 1	Library 37 dilution 1	Library 41 dilution 1	Library 45 dilution 1												
H	Library 1 dilution 2	Library 5 dilution 2	Library 9 dilution 2	Library 13 dilution 2	Library 17 dilution 2	Library 21 dilution 2	Library 25 dilution 2	Library 29 dilution 2	Library 33 dilution 2	Library 37 dilution 2	Library 41 dilution 2	Library 45 dilution 2												

Figure 3. GeneRead Library Quant Array Layout (384-well plate).

- Carefully, tightly seal the plate with Optical Thin-Wall 8-Cap Strips (Formats A and D) or Optical Adhesive Film (Formats C, E, F and G).

**Note:** Users of Bio-Rad and Eppendorf real-time cyclers must ensure that the real-time cycler has been calibrated to use clear, flat optical caps with plates prior to initiating the run.

- Centrifuge for 1 min at 1000 x g at room temperature (15–25°C) to remove bubbles. Visually inspect the plate from underneath to ensure no bubbles are present in the wells.

**Note:** The presence of bubbles in the wells interferes with results.

- Place the plate on ice while setting up the PCR cycling program.

**Note:** The plates containing PCR components mix may be stored at –20°C wrapped in aluminum foil for up to one week if desired.

9. Program the real-time cycler according to Table 3, 4 or 5, depending on the real-time cycler used. If prompted by your cycler software, select "Absolute Quantitation" to begin.

**Table 3. Cycling conditions\* for Applied Biosystems, Bio-Rad<sup>†</sup>, Stratagene and Eppendorf<sup>‡</sup> cyclers**

Cycles	Duration	Temperature	Comments
1	10 min	95°C	HotStarTaq DNA Polymerase is activated by this heating step.
30	15 s	95°C	
	30 s	60°C	
	2 min	72°C	Perform fluorescence data collection.

\* Recommended for the following cyclers: Applied Biosystems models 5700, 7000, 7300, 7500, 7700, 7900HT, StepOnePlus, ViiA 7; Bio-Rad models iCycler, iQ5, MyiQ, MyiQ2, CFX96, CFX384; Stratagene models Mx3000P, Mx3005P, Mx4000P; Eppendorf Mastercycler ep realplex models 2, 2S, 4, 4S.

<sup>†</sup> For Bio-Rad models CFX96 and CFX384: adjust the ramp rate to 1°C/s.

<sup>‡</sup> For Eppendorf Mastercycler ep realplex models 2, 2S, 4, and 4S: for the Silver Thermoblock, adjust the ramp rate to 26%; for the Aluminum Thermoblock, adjust the ramp rate to 35%. Refer to the Instrument Setup Guide at <https://www.qiagen.com/us/shop/assay-technologies/real-time-pcr-and-rt-pcr-reagents/rt2-profiler-pcr-arrays/#resources> for detailed setup instructions.



**Table 4. Cycling conditions for Roche LightCycler 480 II\* Cycler**

Cycles	Duration	Temperature	Comments
1	10 min	95°C	HotStarTaq DNA Polymerase is activated by this heating step.
35	15 s	95°C	
	30 s	60°C	Perform fluorescence data collection.
	2 min	72°C	

\* Recommended for the Roche LightCycler 480 II. If using a Roche LightCycler 480 II, adjust the ramp rate to 1.5°C/s for 96-well plate and 2°C/s for 384-well plate. Refer to the "Instrument Setup Guide" at <https://www.qiagen.com/us/shop/assay-technologies/real-time-pcr-and-rt-pcr-reagents/rt2-profiler-pcr-arrays/#resources> for more information on other required changes to settings for Melt Curve Acquisition.

**Table 5. Cycling conditions for Bio-Rad and Takara cyclers and all other cyclers\***

Cycles	Duration	Temperature	Comments
1	10 min	95°C	HotStarTaq DNA Polymerase is activated by this heating step.
30	15 s	95°C	
	30–40 s	55°C	Perform fluorescence data collection. Different cyclers need different lengths of time to detect the fluorescent signal. Choose the appropriate time for the annealing step (55°C) for your cycler.
	60 s	72°C	

\* Recommended for the following cyclers: Bio-Rad/MJ Research models Chromo4, DNA Engine Opticon, DNA Engine Opticon 2; Takara TP-800; all other cyclers.

10. Place the plate in the real-time cycler. If recommended by the cycler user manual, use a compression pad with plate sealed with optical adhesive film (formats C, E, F and G). Start the run.

## Protocol 2: Real-Time PCR for GeneRead Library Quant Array for Ion Torrent or Illumina (Format R)

### Important points before starting

- Ensure that GeneRead qPCR SYBR Green ROX FAST Mastermix (180850) is used with a Rotor-Gene cycler.
- If dealing with a new library prep without prior measurement of the final concentration, we recommend checking the concentration using BioAnalyzer or QIAxcel and diluting the library to about 1 nM before proceeding with the protocol given below.

### Procedure

#### Preparing sample library dilution

1. Prepare a starting 1:20 dilution of the sample library in dilution buffer (2  $\mu$ l sample library + 38  $\mu$ l dilution buffer).
2. Prepare two working dilutions of sample library, Dilution 1 and Dilution 2, as indicated in Table 6.

**Table 6. Sample library dilution**

Dilution	Library	Dilution buffer
Dilution 1 (1:2000)	2 $\mu$ l of 1:20	198 $\mu$ l
Dilution 2 (1:20000)	5 $\mu$ l of 1:2000	45 $\mu$ l

**Note:** Always use freshly made dilutions.

### Setting up the PCR

3. Briefly centrifuge the GeneRead qPCR SYBR Green ROX FAST Mastermix (10–15 s) to bring the contents to the bottom of the tube.

**Note:** As the GeneRead qPCR SYBR Green ROX FAST Mastermix contains HotStart DNA *Taq* Polymerase, which is active only after heat activation, reactions can be prepared at room temperature (15–25°C).

4. Make PCR mix for each sample dilution (triplicates) or standard / NTC wells as directed in Table 7.


**Table 7. PCR mix for standard or sample libraries\***

Component	Rotor-Disc 100	
	For standard and NTC Wells 1–18 (A1–A12, B1–B6) and R (97–100)	For sample library For dilution 1 and 2 (3 wells for each dilution)
RNase-/DNase-free water	240 µl	27 µl
GeneRead qPCR SYBR Green ROX FAST Mastermix	240 µl	36 µl
Template	–	9 µl
<b>Final volume</b>	<b>480 µl</b>	<b>72 µl</b>

\* For multiple wells, an excess volume is provided to allow for pipetting errors. Perform pipetting steps as precisely as possible to ensure that each well receives the required volume.

5. Slide the plate into the Rotor-Disc 100 Loading Block using the tab at position A1 and the tube guide holes.
6. Add 20 µl per well PCR mix to three wells for each sample dilution (see Figure 4 for layout). Add 20 µl PCR mix (with no template) into each well for well 1–18 (A1 to 12, B1 to B6) and 97–100 (R).  
**Note:** Although wells 97–100 (R) do not contain assays, it is essential to add PCR components mix for optimized balancing.

The layout of figure 4 must be followed as shown for appropriate data analysis.

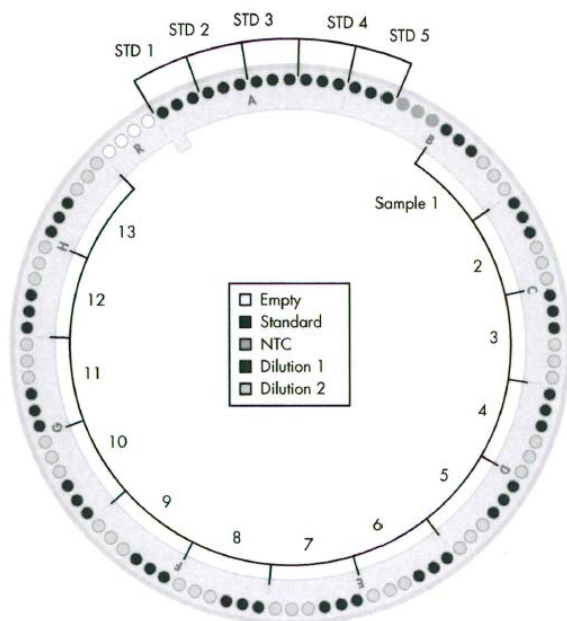

**Figure 4. GeneRead Library Quant Array layout (Format R).**

7. Carefully seal the plate with Rotor-Disc Heat-Sealing Film using the Rotor-Disc Heat Sealer. For detailed instructions, see the *Rotor-Gene Q User Manual*.
8. Program the real-time cycler according to Table 8.

**Note:** For additional help with instrument setup, see our instrument-specific setup instructions and protocol files at:

<https://www.qiagen.com/us/shop/assay-technologies/real-time-pcr-and-rt-pcr-reagents/rt2-profiler-pcr-arrays/#resources>

**Table 8. Cycling conditions for Rotor-Gene cyclers**

Cycles	Duration	Temperature	Comments
1	10 min	95°C	HotStarTaq DNA Polymerase is activated by this heating step
30	15 s	95°C	
	5 s	68°C	
	5 s	65°C	
	60 s	60°C	Perform fluorescence data collection

9. Insert the plate into the Rotor-Disc 100 Rotor and secure with the Rotor-Disc 100 Locking Ring. Start the run. For detailed instructions, see the *Rotor-Gene Q User Manual*.

# Data Analysis

## Obtaining raw threshold cycle ( $C_T$ ) values

After the cycling program has completed, obtain the  $C_T$  values according to the instructions provided by the manufacturer of the real-time PCR instrument. We recommend manually setting the baseline and threshold values according to the following instructions.

Baseline: Use automatic baseline.

Threshold value: Using the Log View of the amplification plots, place the threshold above the background signal but within the lower third of the linear portion of the amplification curves.

## Exporting $C_T$ values

Export and/or copy/paste the  $C_T$  values from the instrument software to a blank Microsoft® Excel spreadsheet according to the manufacturer's instructions for the real-time PCR instrument.

## Excel-based data analysis file

First, download the GeneRead Library Quant Array Excel-based data analysis file that corresponds to the plate format used, which is available at:

**<https://www.qiagen.com/us/shop/sample-technologies/dna-sample-technologies/genomic-dna/generead-library-quant-system/#resources>**

Second, paste the raw  $C_T$ s (for the whole plate, even if you only use part of the plate) into the GeneRead Library Quant Array Excel-based data analysis file.

Enter the required run parameters.

1. Dilution fold
  - a. If you are using the default dilutions (2000 and 20000), no changes are necessary
  - b. If you are not using the default dilutions (2000 and 20000), enter your library dilutions
2. Library size
  - a. If you used the GeneRead DNaseq panels with the GeneRead Library Prep system, do not change the library size in the data analysis file for either Illumina or Ion Torrent.
  - b. If you used a different target enrichment or library construction method, enter the appropriate size of library fragments.

Analyze the automatically generated results by following the directions in the “Instructions” worksheet of the Excel file.

Finally, dilute the original library accordingly to perform emPCR (for Ion Torrent platforms) or dilute the original library to designed range for cluster generation (for Illumina platforms). If you used the GeneRead DNaseq panels with the GeneRead Library Prep system, refer to Appendix F of the all-bead *GeneRead DNaseq Targeted Panels V2 Handbook*.

## Troubleshooting Guide

This troubleshooting guide may be helpful in solving any problems that may arise. For more information, see also the Frequently Asked Questions page at our Technical Support Center: [www.qiagen.com/FAQ/FAQList.aspx](http://www.qiagen.com/FAQ/FAQList.aspx). The scientists in QIAGEN Technical Services are always happy to answer any questions you may have about either the information and protocols in this handbook or sample and assay technologies (for contact information, see back cover or visit [www.qiagen.com](http://www.qiagen.com)).

### Comments and suggestions

---

#### No PCR product is generated

- |  |   |
|--|---|
| a) Assay reagents or master mix are incorrectly stored | Storing assay reagents or master mix at inappropriate temperature for extended periods may cause degradation or reduce their activity and PCR amplification efficiency. |
| b) Incorrect real-time PCR cycling program is used     | Be sure to use the correct cycling program, including 10 minutes at 95°C, to fully activate the hot start enzyme in the GeneRead qPCR SYBR Green Mastermix.             |

#### Signal is detected for no-template controls (NTC)

- |  |  |
|--|--|
| a) Reaction is contaminated by nucleic acids         | Follow recommendations for PCR handling (see page 10 for details).   |
| b) Primer dimers or false amplification is generated | Primer dimers or false amplification may occur in NTCs. As long as the $C_T$ is high (e.g. >29), it is far enough above the $C_T$ range for the actual template that it will have no effect on quantification. |

#### PCR efficiency is out of the range of 90–110%

- |                                   |  |
|-----------------------------------|--|
| a) PCR conditions are not optimal | Ensure that all the reagents are stored properly and reactions are set up correctly. |
| b) PCR program is wrong.          | Pay close attention to the different PCR program for different machines.             |



## Comments and suggestions

---

**$C_T$  of sample library is out of the range of standard curve ( $C_T$  between 6–25, depending on real-time cycler)**

- |                                      |   |
|--------------------------------------|---|
| a) Library concentration is too high | Dilute the sample further to ensure that at least one of the diluted concentrations is within the range of serial dilution standards. |
| b) Library concentration is too low  | Dilute the sample less to ensure at least one of the diluted concentrations is within the range of serial dilution standards.         |

## References

QIAGEN maintains a large, up-to-date online database of scientific publications utilizing QIAGEN products. Comprehensive search options allow you to find the articles you need, either by a simple keyword search or by specifying the application, research area, title, etc.

For a complete list of references, visit QIAGEN's Knowledge Hub online at [www.qiagen.com/knowledgehub](http://www.qiagen.com/knowledgehub) or contact QIAGEN Technical Services or your local distributor.

### Cited references

1. Buehler, B. (2010) Rapid quantification of DNA libraries for next-generation sequencing. *Methods* 50, S15.

## Ordering Information

Product	Contents	Cat. no.
GeneRead Library Quant Array	Two arrays in Formats A, C, D, F, or R; 1 array in Format E or G	180611
Related products		
GeneRead Library Quant Kit	Laboratory-verified forward and reverse primers for 500 x 25 µl reactions (500 µl); Standard (100 µl); Dilution Buffer (30ml); GeneRead qPCR SYBR Green Mastermix (varies)	180612
GeneRead qPCR SYBR Green Mastermix (2)*	For 2 x 96 assays in 96-well plates; suitable for use with real-time cyclers that do not require a reference dye, including: Bio-Rad models CFX96, CFX384, Bio-Rad/MJ Research Chromo4, Bio-Rad/MJ Research Opticon 2; Roche LightCycler 480 II(96-well and 384-well); all other cyclers	180820
GeneRead qPCR SYBR Green ROX Mastermix (2)*	For 2 x 96 assays in 96-well plates; suitable for use with the following real-time cyclers: Applied Biosystems models 5700, 7000, 7300, 7500 [Standard and Fast], 7700, 7900HT 96-well block [Standard and Fast] and 384-well block, StepOnePlus; Eppendorf Mastercycler ep realplex models 2, 2S, 4, 4S; Stratagene models Mx3000P, Mx3005P, Mx4000	180840
GeneRead qPCR SYBR Green Fluor Mastermix (2)*	For 2 x 96 assays in 96-well plates; suitable for use with the following real-time cyclers: Bio-Rad models iCycler, iQ5, MyiQ, MyiQ2	180830

\* Larger kit sizes available; please inquire.

Product	Contents	Cat. no.
GeneRead qPCR SYBR Green ROX FAST Mastermix (2)*	For 2 x 96 assays, suitable for use with Rotor-Gene Q and Rotor-Gene 6000	180850

For up-to-date licensing information and product-specific disclaimers, see the respective QIAGEN kit handbook or user manual. QIAGEN kit handbooks and user manuals are available at [www.qiagen.com](http://www.qiagen.com) or can be requested from QIAGEN Technical Services or your local distributor.

**Notes**

**Notes**



Trademarks: QIAGEN®, GeneRead™ HotStarTaq®, Rotor-Gene®, QIAGEN Group); Agilent®, BIAanalyzer® (Agilent Technologies, Inc.); Applied Biosystems®, Ion Torrent™, StepOnePlus™, ROX™, SYBR® (Life Technologies Corporation); Bio-Rad®, iCycler®, Chromo4™, CFX96™, DNA Engine Opticon®, CFX384™, iQ™, MyiQ™, MyiQ2™ (Bio-Rad Laboratories, Inc.); Eppendorf®, Mastercycler® (Eppendorf AG); Illumina® (Illumina, Inc.); Roche®, LightCycler® (Roche Group); SmartCycler® (Cepheid); Stratagene®, Mx3005P®, Mx3000P®, Mx4000® (Stratagene); Excel®, Microsoft® (Microsoft Corporation). Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

Purchaser agrees that use of this product and data therefrom is limited solely to the purchaser and for only the purchaser's own internal molecular biology research applications ("Permitted Use"), and shall not be re-sold or used for any other purposes (all of which are expressly prohibited), including without limitation diagnostic purposes, uses that could require regulatory approval for diagnostics from an agency of any government or regulatory entity anywhere in the world, diagnosis, prevention, or treatment of disease, and the right to perform commercial services of any kind, including without limitation, reporting the results of purchaser's activities, including without limitation, for a fee or other commercial consideration. Except for the Permitted Use, no rights, titles, or interests in or to any tangible or intangible property rights are conveyed or shall be deemed conveyed by implication, estoppel or otherwise. The performance characteristics of the product other than for the Permitted Use are unknown.

#### Limited License Agreement

Use of this product signifies the agreement of any purchaser or user of the GeneRead Library Quant Array Handbook to the following terms:

1. The Ion Torrent Library Quantification Array may be used solely in accordance with the *GeneRead Library Quant Array Handbook* and for use with components contained in the Kit only. QIAGEN grants no license under any of its intellectual property to use or incorporate the enclosed components of this Kit with any components not included within this Kit except as described in the *GeneRead Library Quant Array Handbook* and additional protocols available at [www.qiagen.com](http://www.qiagen.com).
2. Other than expressly stated licenses, QIAGEN makes no warranty that this Kit and/or its use(s) do not infringe the rights of third-parties.
3. This Kit and its components are licensed for one-time use and may not be reused, refurbished, or resold.
4. QIAGEN specifically disclaims any other licenses, expressed or implied other than those expressly stated.
5. The purchaser and user of the Kit agree not to take or permit anyone else to take any steps that could lead to or facilitate any acts prohibited above. QIAGEN may enforce the prohibitions of this Limited License Agreement in any Court, and shall recover all its investigative and Court costs, including attorney fees, in any action to enforce this Limited License Agreement or any of its intellectual property rights relating to the Kit and/or its components.

For updated license terms, see [www.qiagen.com](http://www.qiagen.com).

This product is compatible for use in real-time PCR methodologies, including 5' nuclease processes and dsDNA-binding dye processes claimed in patents owned by Roche Molecular Systems, Inc., F. Hoffmann-La Roche Ltd, and Applied Biosystems Corporation. No right to perform any of these patented processes is conveyed, expressly, by implication or by estoppel to the purchaser by the purchase of this product. A license to use these patented processes for the purchaser's own internal research accompanies the purchase of certain reagents of Applied Biosystems and other licensed suppliers, or is available from Applied Biosystems. Further information on purchasing licenses may be obtained by contacting the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.

© 2015 QIAGEN, all rights reserved.

---

[www.qiagen.com](http://www.qiagen.com)

**Australia** ■ Orders 1-800-243-800 ■ Fax 03-9840-9888 ■ Technical 1-800-243-066

**Austria** ■ Orders 0800-28-10-10 ■ Fax 0800-28-10-19 ■ Technical 0800-28-10-11

**Belgium** ■ Orders 0800-79612 ■ Fax 0800-79611 ■ Technical 0800-79556

**Brazil** ■ Orders 0800-557779 ■ Fax 55-11-5079-4001 ■ Technical 0800-557779

**Canada** ■ Orders 800-572-9613 ■ Fax 800-713-5951 ■ Technical 800-DNA-PREP (800-362-7737)

**China** ■ Orders 86-21-3865-3865 ■ Fax 86-21-3865-3965 ■ Technical 800-988-0325

**Denmark** ■ Orders 80-885945 ■ Fax 80-885944 ■ Technical 80-885942

**Finland** ■ Orders 0800-914416 ■ Fax 0800-914415 ■ Technical 0800-914413

**France** ■ Orders 01-60-920-926 ■ Fax 01-60-920-925 ■ Technical 01-60-920-930 ■ Offers 01-60-920-928

**Germany** ■ Orders 02103-29-12000 ■ Fax 02103-29-22000 ■ Technical 02103-29-12400

**Hong Kong** ■ Orders 800 933 965 ■ Fax 800 930 439 ■ Technical 800 930 425

**Ireland** ■ Orders 1800 555 049 ■ Fax 1800 555 048 ■ Technical 1800 555 061

**Italy** ■ Orders 800-789-544 ■ Fax 02-334304-826 ■ Technical 800-787980

**Japan** ■ Telephone 03-6890-7300 ■ Fax 03-5547-0818 ■ Technical 03-6890-7300

**Korea (South)** ■ Orders 080-000-7146 ■ Fax 02-2626-5703 ■ Technical 080-000-7145

**Luxembourg** ■ Orders 8002-2076 ■ Fax 8002-2073 ■ Technical 8002-2067

**Mexico** ■ Orders 01-800-7742-639 ■ Fax 01-800-1122-330 ■ Technical 01-800-7742-436

**The Netherlands** ■ Orders 0800-0229592 ■ Fax 0800-0229593 ■ Technical 0800-0229602

**Norway** ■ Orders 800-18859 ■ Fax 800-18817 ■ Technical 800-18712

**Singapore** ■ Orders 1800-742-4362 ■ Fax 65-6854-8184 ■ Technical 1800-742-4368

**Spain** ■ Orders 91-630-7050 ■ Fax 91-630-5145 ■ Technical 91-630-7050

**Sweden** ■ Orders 020-790282 ■ Fax 020-790582 ■ Technical 020-798328

**Switzerland** ■ Orders 055-254-22-11 ■ Fax 055-254-22-13 ■ Technical 055-254-22-12

**UK** ■ Orders 01293-422-911 ■ Fax 01293-422-922 ■ Technical 01293-422-999

**USA** ■ Orders 800-426-8157 ■ Fax 800-718-2056 ■ Technical 800-DNA-PREP (800-362-7737)