



October 2024

# QIAcuity<sup>®</sup> Software Suite Backup and Restore Scripts User Guide

For QIAcuity Software Suite version 3.0

# Table of Contents

1. Introduction .....	3
2. Preconditions .....	4
3. Start the Backup Manually .....	5
4. Setup an Automated Backup with Microsoft® Task Scheduler .....	6
5. Restore a Backup .....	13
6. Troubleshooting .....	15
6.1. Inconsistency of regional settings during restoring from backup .....	16
Document Revision History .....	17

# 1. Introduction

The QIAcuity Software Suite Backup and Restore Scripts for Software Suite 3.0 are a standalone solution to backup all relevant user data of the QIAGEN® Software Suite for disaster recovery and to restore the data on a new or existing QIAGEN Software Suite installation. The scripts are designed for QIAcuity Software Suite version 3.0 only. The backup can be conducted manually or automated by using of the Windows® task scheduler. In detail, the following data are backed up:

1. User Management Database
2. Audit Trail Database and Tracker setting
3. Volume Precision Factor (VPF)
4. Plates
5. Plate Templates
6. Reports
7. Archive configurations (archive location, automatic archiving setting)
8. Instrument registration list

The plate archive is not part of the backup.

**Note:** For backup and restore of Software Suite versions lower than 3.0, please refer to the dedicated QIAcuity Software Suite Backup and Restore Scripts for versions 2.5.0.0, 2.5.0.1, and 2.5.1.0, or for version 2.0, 2.1.7, 2.1.8, and 2.2.

## 2. Preconditions

Windows administrator permission is needed to setup and perform an automated backup and for manual backup and restore.

Visit [www.qiagen.com](http://www.qiagen.com) and go to the Latest Software Version section of the QIAcuity product page to download the backup and restore package for QIAcuity Software Suite version 3.0.

Unzip the downloaded package to a location of your choice. Please note that **suite\_backup.exe**, **suite\_backup.ps1**, and **suite\_rollback.exe** must be in the same directory.

To access network folders for backup and restore, the UNC (Universal Naming Convention, `\\servername\path`) path format has to be used.

For example, if you have a mapped network drive on the source PC of your backup, you can identify the UNC path by using the command "net use" in a PowerShell® or CMD command prompt. If you do not know how to do this, please contact your local IT administrator.

**Important:** To ensure data integrity of the databases during backup and restore, the QIAcuity Software Suite is terminated for backup and restore. No plate runs, analyses, archiving, or any other functionality will be available during backup or restore. Please make sure that the backup and restore run only if no plate run and analysis or automatic archiving is in progress; otherwise, the data will be lost. If automatic archiving is activated, it will start at 1 AM Suite PC time and will be skipped if the QIAcuity Software Suite is still terminated for the backup or restore at that time.

**Note:** The time required for a backup depends on various factors such as the size of the databases, the speed of the interface of the backup device (USB, network, etc.), and the utilization of the PC where the backup or restore script was started.

You can estimate the duration of the backup by using the following equation:

$$\text{Estimated Backup duration} = \frac{\text{Size of Suite Directory} + \text{Size of PostgreSQL Data Directory}}{\text{System performance}}$$

**Suite Directory** = C:\Program Files (x86)\QIAcuity Software Suite

**PostgreSQL Data Directory** = C:\ProgramData\QIAcuity\SuiteDbData


As an example, assuming an overall available system performance of 70 Mbytes/s and a Suite Directory + PostgreSQL Data Directory size of 100 GByte, the backup takes approximately 24 minutes.

### 3. Start the Backup Manually

Open the backup and restore directory with the file explorer and run **suite\_backup.exe** as administrator.

The backup script starts with opening a command shell window showing the output of the script that is also logged into the **data\_backup.log** file. The standard setting for the Suite directory and backup file target drive and directory (**C:\**) are used.

If you want to change the target drive or directory, you can run the script from the command line providing the target as a parameter. For example, the target is a local drive and directory **D:\Backup**, then the command line is "suite\_backup D:\Backup".

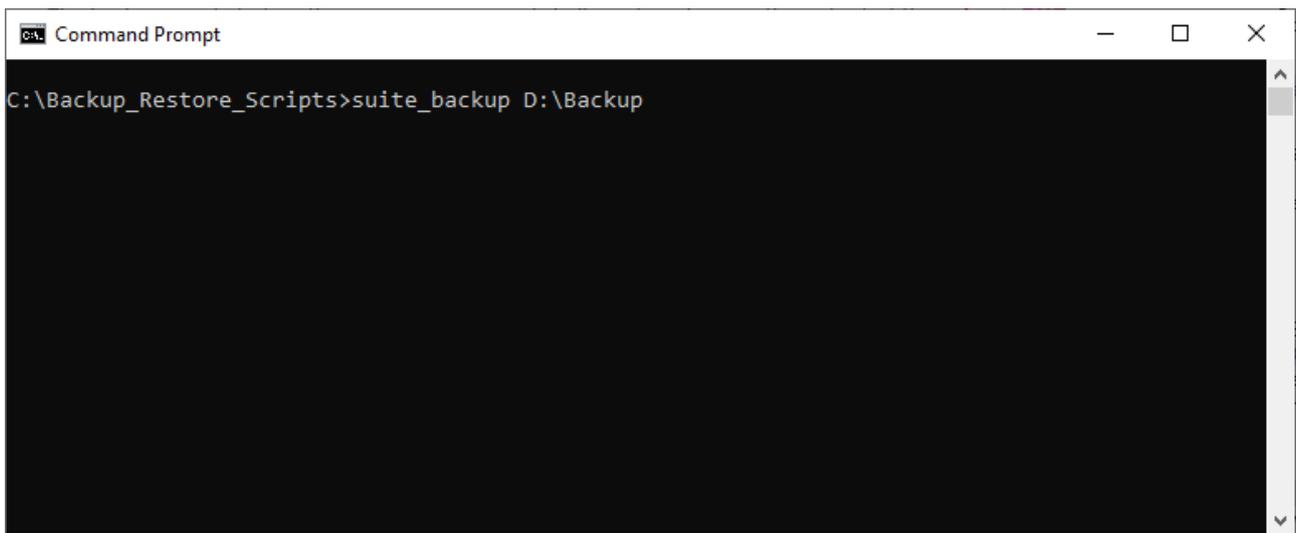

If the target is a network share, the UNC convention for the target has to be used; the command line is "suite\_backup\\backupserver\_name\backupdirectory".


**Note:** When typing in the command lines, enter the exact text described above, including all spaces in between the command.

In each case, the backup is being stored in a subdirectory **suite\_backup\_%timestamp%** of your backup directory.


**Important:** Following the best practices for backing up data, ensure that each backup is restorable before archiving it. If you do not know how to do this, please contact your local IT department.

## 4. Setup an Automated Backup with Microsoft® Task Scheduler


1. Search for the windows task scheduler application and run it.

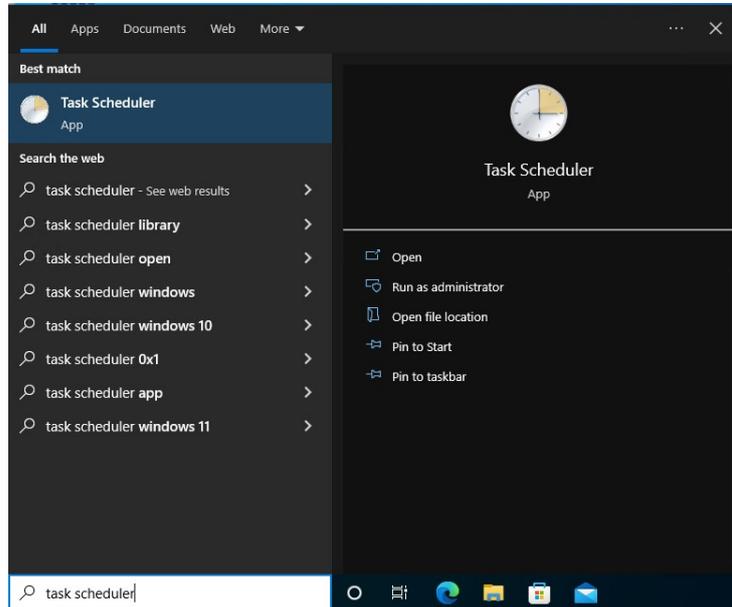


2. Select **Create Task...** action from right menu (highlighted).

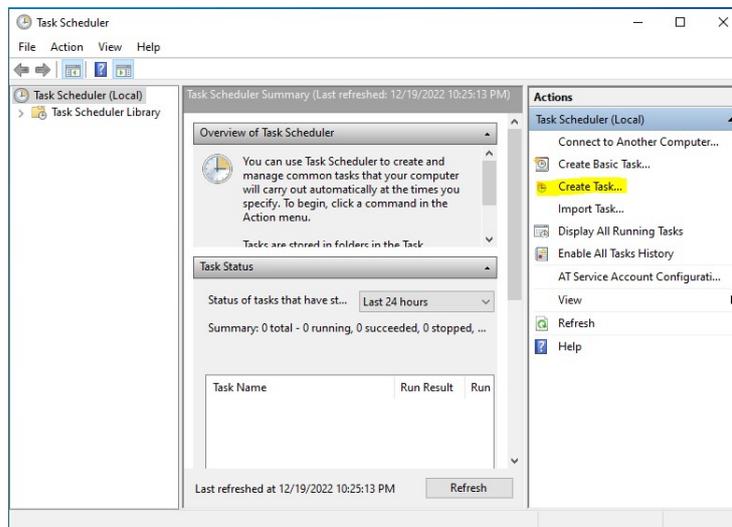


3. In the **General** tab, choose a name for the task.

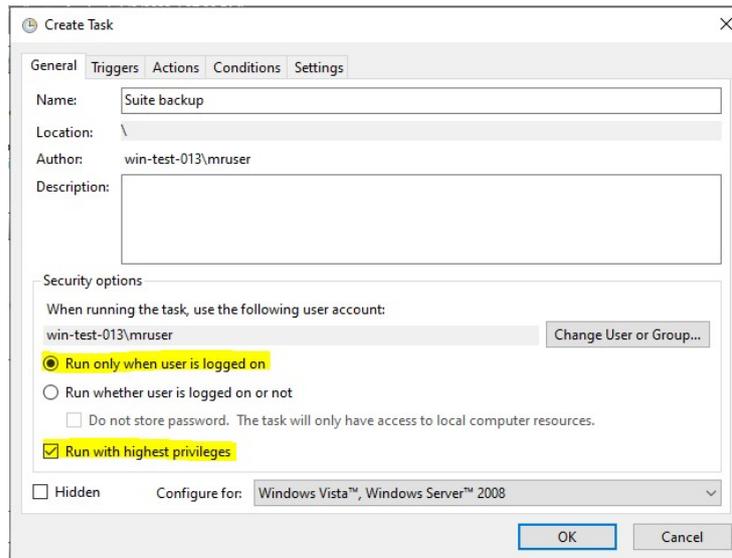

4. In **Security options** of the Create Task window, you will have to choose whether you want to see error messages (in case of an error) like this or not:

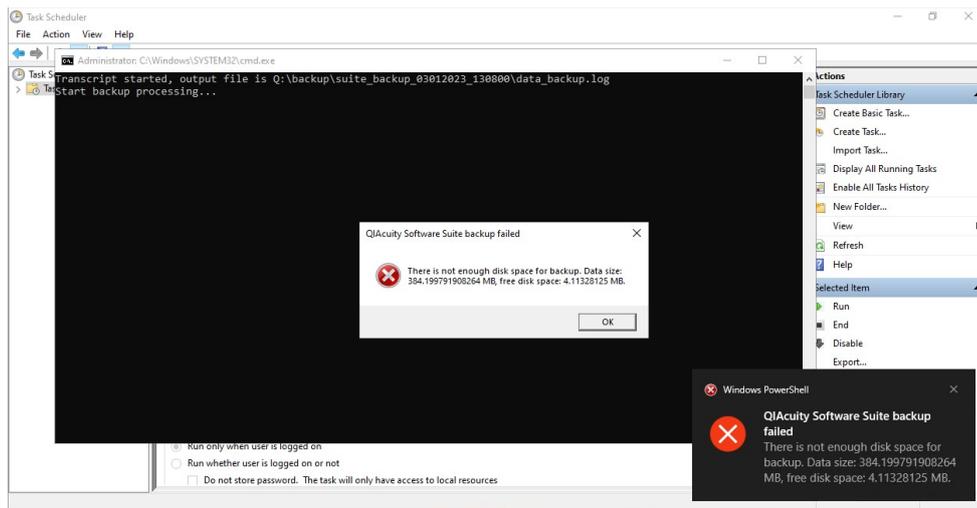


- a. If you want to run the backup only when a user is logged on, to be informed of errors in that manner, pick the options as highlighted in the screenshot below.

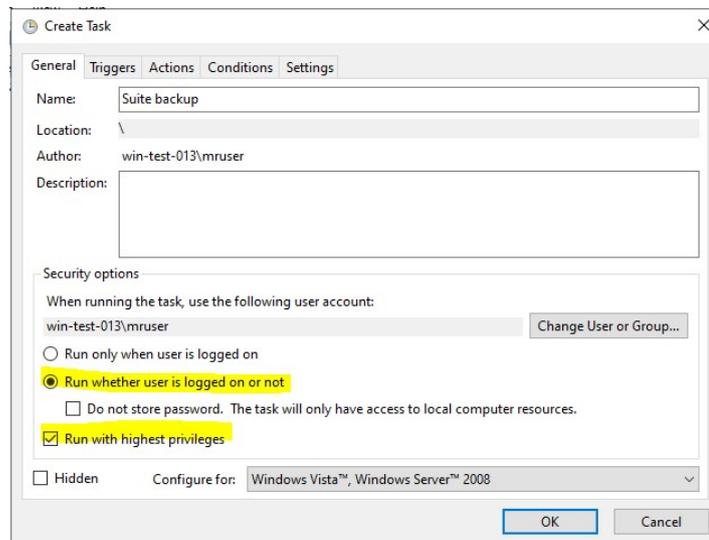


The backup process will be performed only when you are logged in, and you will see the error messages as presented in the screenshot above. In addition, any error will be logged in the backup log file.


- b. If you want to run the backup process in the background without a user logged in, pick the options as highlighted in the screenshot below. With this option, any backup errors will be logged in the backup log file only.

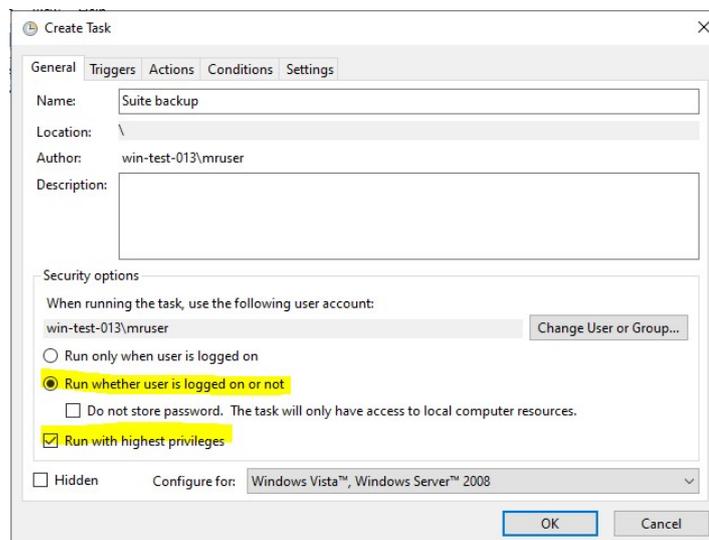


5. In the **Triggers** tab, click the **New...** button.

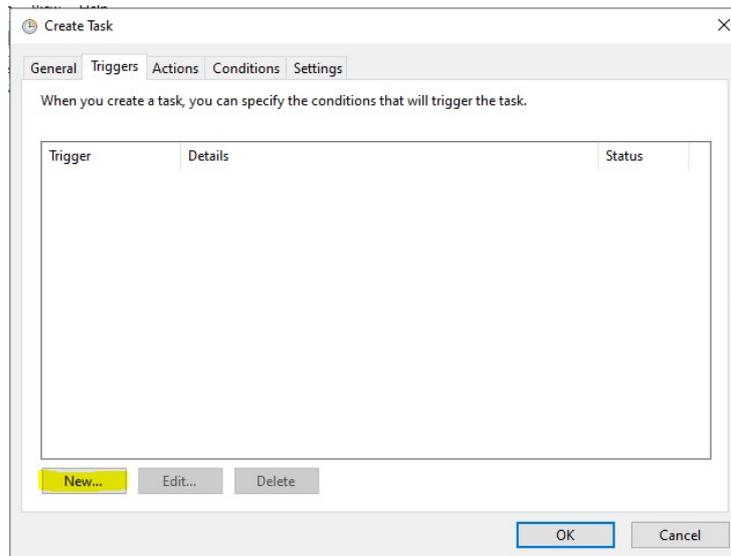

The New Trigger configuration window allows to define when the backup action should be performed.

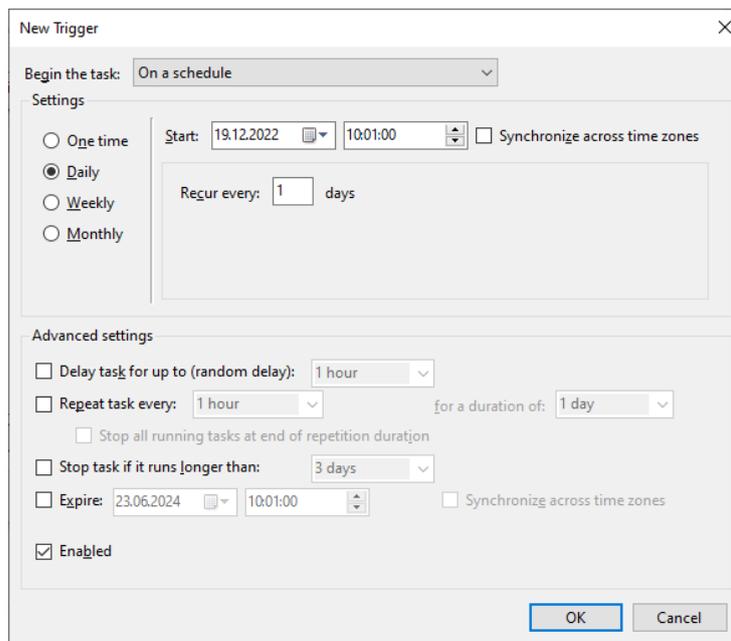


6. Go to the **Actions** tab, and click the **New...** button:

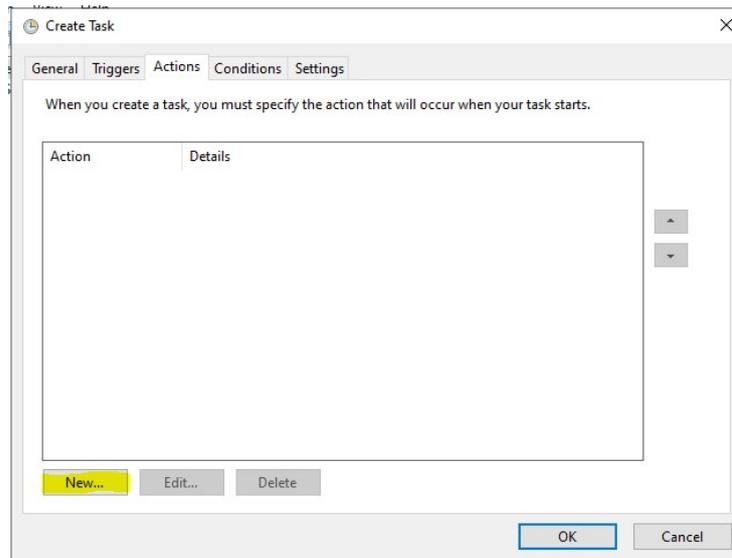


7. In the New Action window, you have to provide the path to the backup script on your disc ("Program/script" text area) and your desired, existing directory to store the automatic backups ("Add arguments (optional)" text area).


If the target is a network share, the UNC convention for the argument has to be used; for example, `\\backupserver_name\backupdirectory`.

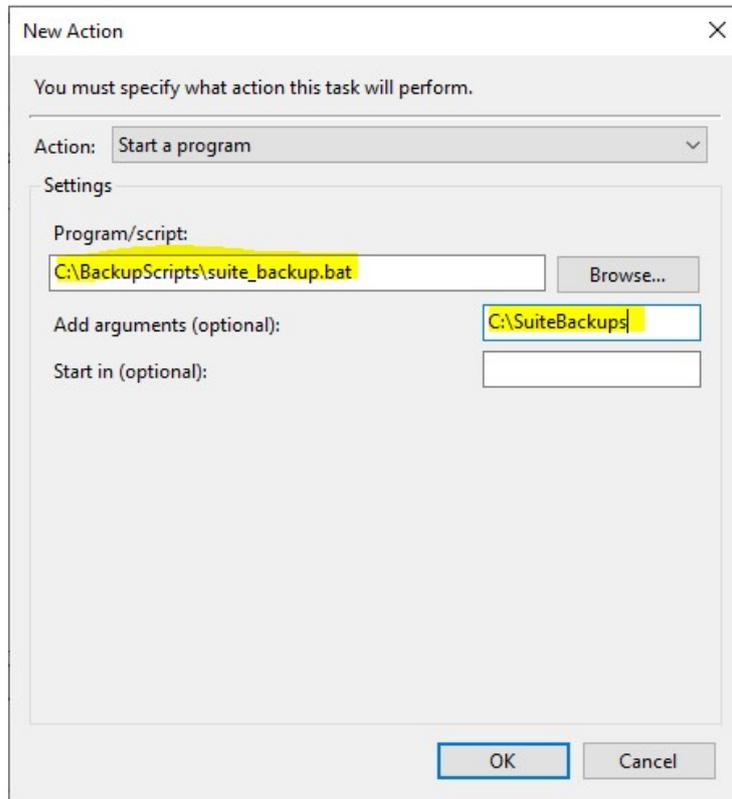


8. Configuration in the **Conditions** and **Settings** tabs is up to your needs. For the most common case, nothing needs to be adjusted from the default.

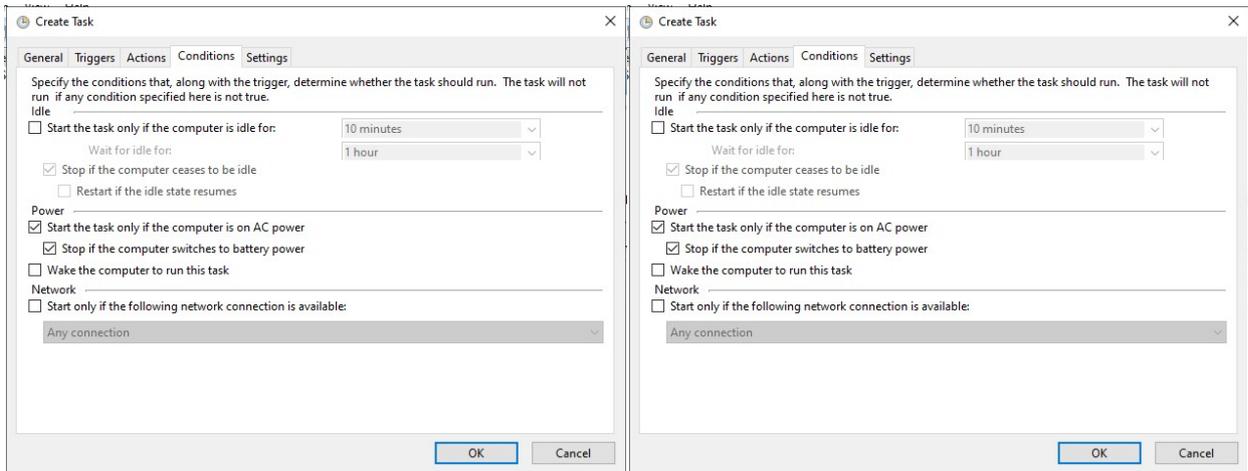

9. With **Triggers** and **Actions** tab filled, click the **OK** button. In case the configuration to run the backup in the background was selected, Windows will ask you for the password of the user account running the task. Provide it, and click **OK**.

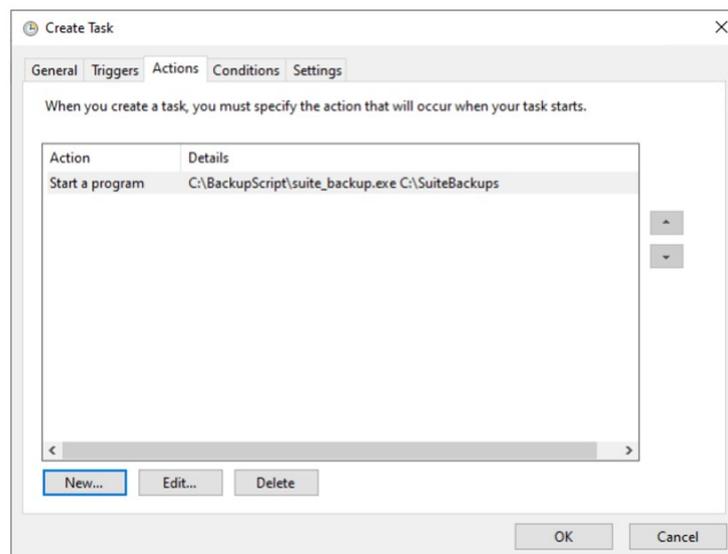


After that step, you can see your task in the Task Scheduler Library view.

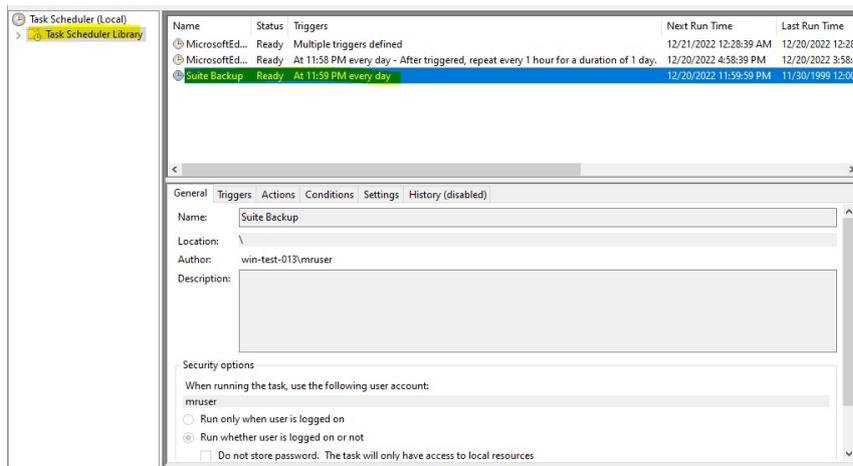

The windows task scheduler will from now on automatically start the Suite Backup Task and create the Suite backups as configured in the task.

Inside each timestamped backup folder (**suite\_backup\_%timestamp%**), the backed-up files along with a backup log file and a rollback script are stored.

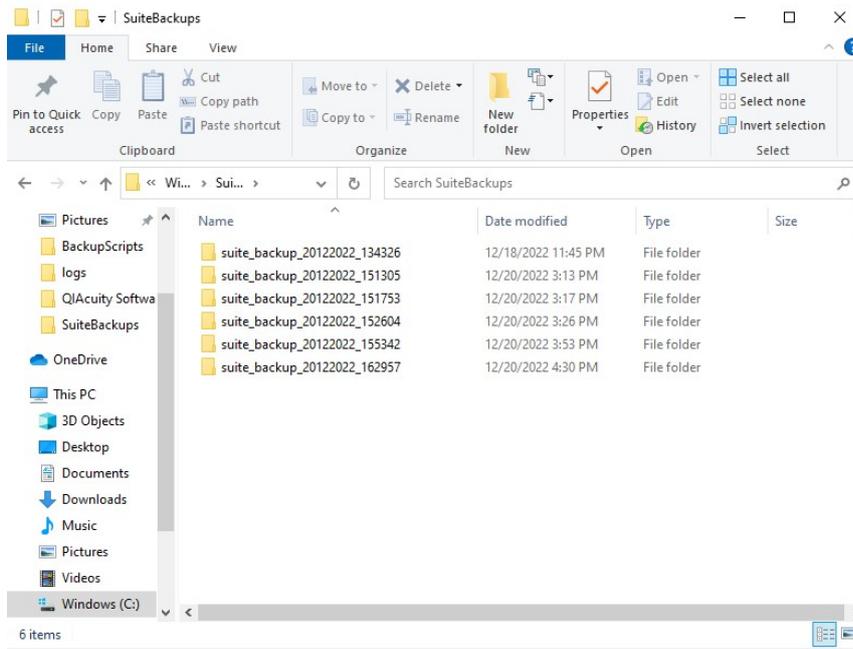


## 5. Restore a Backup


**Important:** In case of a disaster recovery, install the QIAcuity Software Suite first. Please make sure that the software version is identical to the version of the backup.

**Important:** The rollback script will restore the data of the backup to the QIAcuity Software Suite Installation directories and overwrite existing data.


To restore a backup, run **suite\_rollback.exe** as administrator inside the desired backup folder and follow the instructions.

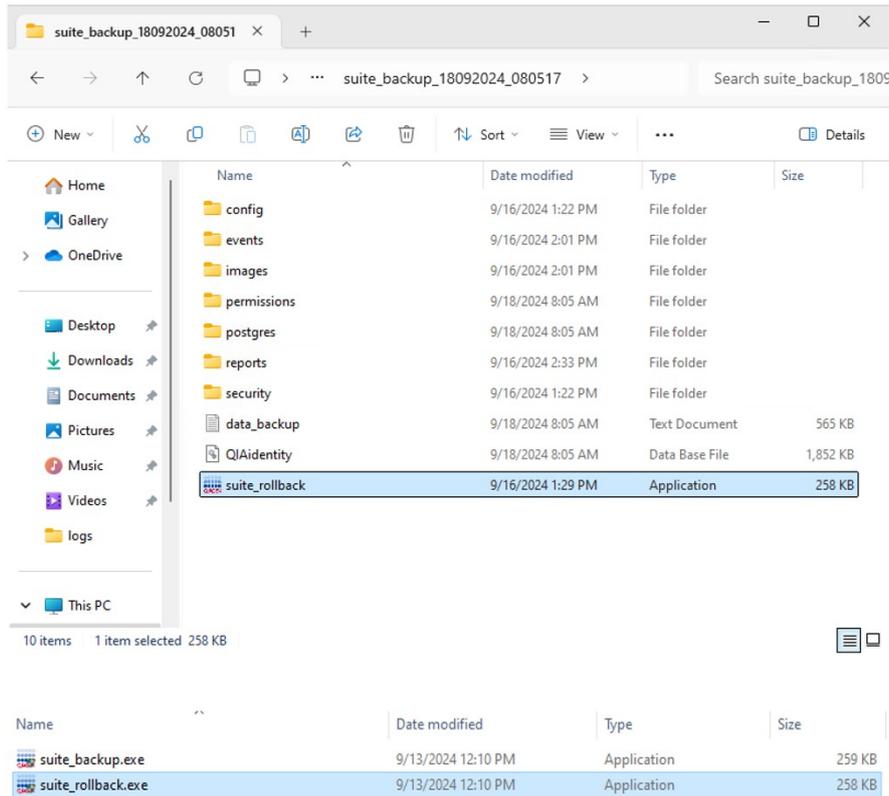

Click **Yes** on below pop-up window.

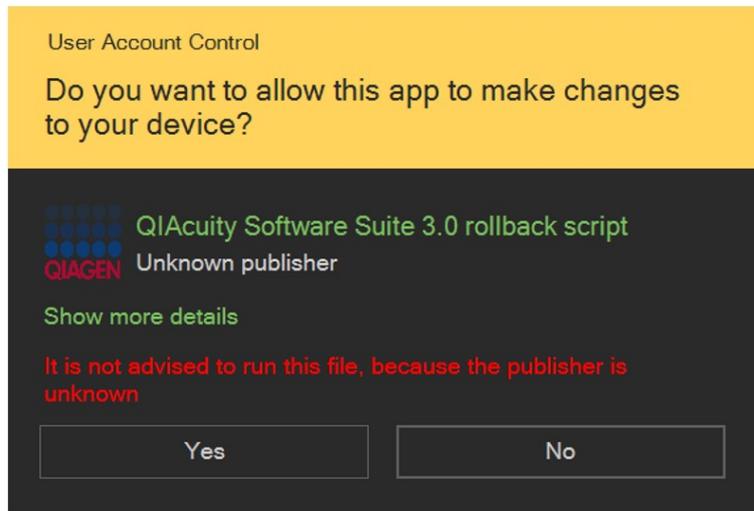

Confirm that the backup data can be restored to the target location by entering "rollback" into the textbox. This will overwrite the data at the target location.

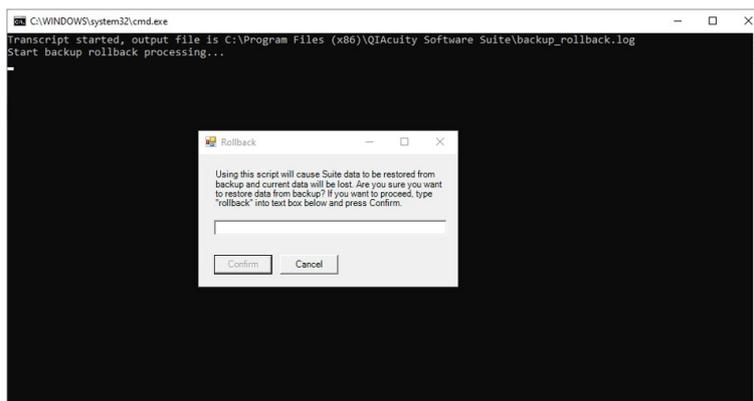


After the rollback has finished, the command window closes and the QIAcuity Software Suite can be opened again.

# 6. Troubleshooting

The backup and restore scripts are logging the progress and results in log files. For the backup script, the logged information is located in ...\**suite\_backup\_%timestamp%\ data\_backup.log**.


For the restore script, the logged information is located in **C:\Program Files (x86)\QIAcuity Software Suite\backup\_rollback.log**.

In case of an error, like incomplete backup, empty backup directory, or incomplete restore, please check the log files for errors.

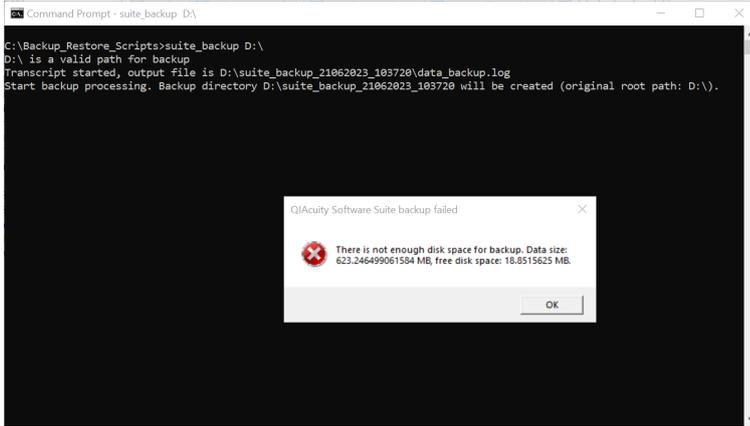Error	Solution
<pre>PS&gt;TerminatingError(): "service QIAcuitySuite cannot be stopped" &gt;&gt; TerminatingError(): "service QIAcuitySuite cannot be stopped" Backup failed: service QIAcuitySuite cannot be stopped</pre>	Run the backup script with administrator rights.
<p>The command line window closes itself right after the script gets started and the backup directory is not created on the network drive, no log is present</p>	Check that the user performing the backup has permission to write on the network drive.
<p>Error message pop up during manual backup about not enough disk space. The command line remains open.</p>	Check the disk space at the backup target location
	
<p>For automated backup, the error message "PS&gt;TerminatingError(): "There is not enough disk space for backup." is present in the data_backup.log file, resulting in no files in the backup directory</p>	
<p>Backup restore made by Windows user different than the one who installed Suite causes internal error while trying to log in to the application.</p>	Restore the backup one more time.

## 6.1. Inconsistency of regional settings during restoring from backup

In case the following pop-up window appears during backup restoration, change the regional setting of the operating system to meet the regional settings that was set during installation of backed-up QIAcuity Software Suite, as stated in the warning.

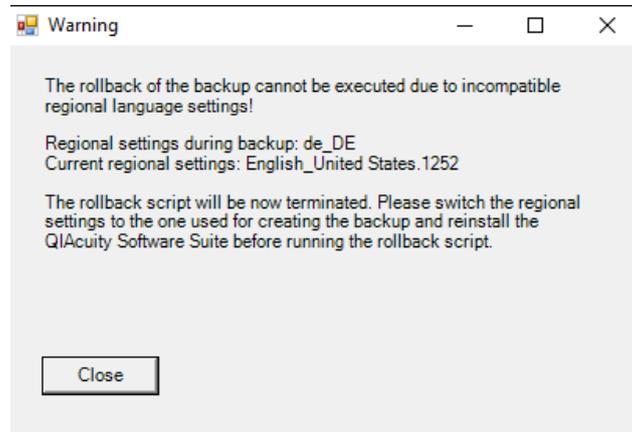

# Document Revision History

Revision	Description
October 2024	Initial release for Software Suite version 3.0.

For up-to-date licensing information and product-specific disclaimers, see the respective QIAGEN kit handbook or user manual. QIAGEN kit handbooks and user manuals are available at [www.qiagen.com](http://www.qiagen.com) or can be requested from QIAGEN Technical Services or your local distributor.

Trademarks: QIAGEN®, Sample to Insight®, QIAcuity® (QIAGEN Group); Microsoft®, PowerShell®, Windows® (Microsoft Corporation). Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

HB-3636-001 10/2024 © 2024 QIAGEN, all rights reserved.

Page left intentionally blank.

