

QIAamp[®] DNA Micro Handbook

For purification of genomic DNA from

small volumes of blood

dried blood spots

urine

tissues

laser-microdissected tissues

and for cleanup of genomic DNA



QIAGEN Sample and Assay Technologies

QIAGEN is the leading provider of innovative sample and assay technologies, enabling the isolation and detection of contents of any biological sample. Our advanced, high-quality products and services ensure success from sample to result.

QIAGEN sets standards in:

- Purification of DNA, RNA, and proteins
- Nucleic acid and protein assays
- microRNA research and RNAi
- Automation of sample and assay technologies

Our mission is to enable you to achieve outstanding success and breakthroughs. For more information, visit www.qiagen.com.



QIAGEN is a member of the Forest Stewardship Council (FSC). For the production of printed materials, including handbooks, QIAGEN has a policy to select suppliers that comply with FSC standards for printing processes and well-managed forests.

Contents

Kit Contents	4
Storage	4
Product Use Limitations	5
Product Warranty and Satisfaction Guarantee	5
Technical Assistance	5
Safety Information	6
Quality Control	7
Introduction	7
Principle and procedure	7
Automated purification	10
Equipment and Reagents to Be Supplied by User	11
Important Notes	12
Eluting pure DNA	12
Yield and size of DNA	12
Carrier RNA	13
Handling of QIAamp MinElute columns	13
Preparation of buffers	14
Protocols	
■ Isolation of Genomic DNA from Small Volumes of Blood	16
■ Isolation of Genomic DNA from Dried Blood Spots	19
■ Isolation of Genomic DNA from Urine	22
■ Isolation of Genomic DNA from Tissues	25
■ Isolation of Genomic DNA from Laser-Microdissected Tissues	28
■ Cleanup of Genomic DNA	31
Troubleshooting Guide	33
Appendix: Working with DNA	35
Ordering Information	36

Kit Contents

QIAamp® DNA Micro Kit	(50)
Catalog no.	56304
Number of preps	50
QIAamp MinElute™ Columns	50
Collection Tubes (2 ml)	200
Buffer ATL	50 ml
Buffer AL*	33 ml
Buffer AW1* (concentrate)	19 ml
Buffer AW2† (concentrate)	13 ml
Buffer AE	12 ml
Carrier RNA (red cap)	310 µg
Proteinase K	1.25 ml
Handbook	1

* Contains a guanidine salt. Not compatible with disinfectants containing bleach. See page 6 for safety information.

† Contains sodium azide as a preservative.

Storage

QIAamp MinElute columns should be stored at 2–8°C upon arrival. However, short-term storage (up to 4 weeks) at room temperature (15–25°C) does not affect their performance.

All buffers can be stored at room temperature (15–25°C) for up to 1 year without showing any reduction in performance.

QIAamp DNA Micro Kits contain a novel, ready-to-use proteinase K solution, which is supplied in a specially formulated storage buffer. The proteinase K is stable for up to 1 year after delivery when stored at room temperature (15–25°C). For storage longer than one year or if ambient temperatures often exceed 25°C, we suggest storing proteinase K at 2–8°C.

Product Use Limitations

The QIAamp DNA Micro Kit is intended for molecular biology applications. This product is not intended for the diagnosis, prevention, or treatment of a disease.

All due care and attention should be exercised in the handling of the products. We recommend all users of QIAGEN products to adhere to the NIH guidelines that have been developed for recombinant DNA experiments, or to other applicable guidelines.

Product Warranty and Satisfaction Guarantee

QIAGEN guarantees the performance of all products in the manner described in our product literature. The purchaser must determine the suitability of the product for its particular use. Should any product fail to perform satisfactorily due to any reason other than misuse, QIAGEN will replace it free of charge or refund the purchase price. We reserve the right to change, alter, or modify any product to enhance its performance and design. If a QIAGEN product does not meet your expectations, simply call your local Technical Service Department or distributor. We will credit your account or exchange the product — as you wish. Separate conditions apply to QIAGEN scientific instruments, service products, and to products shipped on dry ice. Please inquire for more information.

A copy of QIAGEN terms and conditions can be obtained on request, and is also provided on the back of our invoices. If you have questions about product specifications or performance, please call QIAGEN Technical Services or your local distributor (see back cover or visit www.qiagen.com).

Technical Assistance

At QIAGEN we pride ourselves on the quality and availability of our technical support. Our Technical Service Departments are staffed by experienced scientists with extensive practical and theoretical expertise in sample and assay technologies and the use of QIAGEN products. If you have any questions or experience any difficulties regarding the QIAamp DNA Micro Kit or QIAGEN products in general, please do not hesitate to contact us.

QIAGEN customers are a major source of information regarding advanced or specialized uses of our products. This information is helpful to other scientists as well as to the researchers at QIAGEN. We therefore encourage you to contact us if you have any suggestions about product performance or new applications and techniques.

For technical assistance and more information, please see our Technical Support Center at www.qiagen.com/Support or call one of the QIAGEN Technical Service Departments or local distributors (see back cover or visit www.qiagen.com).

Safety Information

When working with chemicals, always wear a suitable lab coat, disposable gloves, and protective goggles. For more information, please consult the appropriate material safety data sheets (MSDSs). These are available online in convenient and compact PDF format at www.qiagen.com/Support/MSDS.aspx where you can find, view, and print the MSDS for each QIAGEN kit and kit component.



CAUTION: DO NOT add bleach or acidic solutions directly to the sample preparation waste.

Buffers AL and AW1 contain guanidine hydrochloride, which can form highly reactive compounds when combined with bleach.

If liquid containing these buffers is spilt, clean with suitable laboratory detergent and water. If the spilt liquid contains potentially infectious agents, clean the affected area first with laboratory detergent and water, and then with 1% (v/v) sodium hypochlorite. The following risk and safety phrases apply to components of the QIAamp DNA Micro Kit:

Buffer AL and Buffer AW1

Contains guanidine hydrochloride: harmful, irritant. Risk and safety phrases:* R22-36/38, S13-26-36-46

Proteinase K

Contains proteinase K: sensitizer, irritant. Risk and safety phrases:* R36/37/38-42/43, S23-24-26-36/37

24-hour emergency information

Emergency medical information in English, French, and German can be obtained 24 hours a day from:

Poison Information Center Mainz, Germany

Tel: +49-6131-19240

* R22: Harmful if swallowed; R36/38: Irritating to eyes and skin; R36/37/38: Irritating to eyes, respiratory system, and skin; R42/43: May cause sensitization by inhalation and skin contact; S13: Keep away from food, drink, and animal feeding stuffs; S23: Do not breathe spray; S24: Avoid contact with skin; S26: In case of contact with eyes, rinse immediately with plenty of water and seek medical advice; S36: Wear suitable protective clothing; S36/37: Wear suitable protective clothing and gloves; S46: If swallowed, seek medical advice immediately and show this container or label.

Quality Control

In accordance with QIAGEN's ISO-certified Quality Management System, each lot of QIAamp DNA Micro Kits is tested against predetermined specifications to ensure consistent product quality.

Introduction

The QIAamp DNA Micro Kit uses well-established technology for purification of genomic and mitochondrial DNA from small sample volumes or sizes. The kit combines the selective binding properties of a silica-based membrane with flexible elution volumes of between 20 and 100 μ l. The procedure is suitable for a wide range of sample materials such as small volumes of blood, blood cards and small tissue samples.

The procedure is designed to ensure that there is no sample-to-sample cross-contamination and allows safe handling of potentially infectious samples. After sample lysis, the simple QIAamp DNA Micro procedure, which is highly suited for simultaneous processing of multiple samples, yields pure DNA in less than 30 minutes.

DNA is eluted in Buffer AE or water and is immediately ready for use in amplification reactions or for storage at -20°C . Purified DNA is free of proteins, nucleases, and other impurities.

Principle and procedure

The QIAamp DNA Micro procedure consists of 4 steps (see flowchart, page 9):

- Lyse — the sample is lysed
- Bind — the DNA in the lysate binds to the membrane of the QIAamp MinElute column
- Wash — the membrane is washed
- Elute — DNA is eluted from the membrane

Sample volumes

The QIAamp DNA Micro Kit is designed for use with small amounts of sample. Protocols are provided for processing:

- Whole blood samples (1–100 μ l)
- Blood card punches (3 mm diameter)
- Urine (up to 10 ml)
- Small amounts of tissue (<10 mg)
- Laser-microdissected tissues

In addition, the QIAamp DNA Micro Kit can be used to clean up genomic DNA samples of up to 100 µl containing up to 10 µg DNA.

Lysing samples

Samples are lysed under highly denaturing conditions at elevated temperatures in the presence of proteinase K and Buffer ATL.

Binding DNA to the QIAamp MinElute column membrane

To allow optimal binding of DNA to the membrane, Buffer AL (and in most protocols, ethanol) is added to the lysate. Lysates are transferred onto a QIAamp MinElute column, where DNA is adsorbed onto the silica-gel membrane as the lysate is drawn through by centrifugation. Salt and pH conditions ensure that proteins and other contaminants, which can inhibit PCR and other downstream enzymatic reactions, are not retained on the QIAamp MinElute column membrane. DNA yield depends on the volume or size and on the quality of the starting sample.

For cleanup of genomic DNA, the conditions for binding to the silica-gel membrane are adjusted by adding Buffers AW1 and AW2.

QIAamp MinElute columns fit into most standard microcentrifuge tubes. Due to the volume of flow-through, 2 ml collection tubes (provided) are required to support the QIAamp MinElute column during sample loading and wash steps.


Washing away residual contaminants

While nucleic acids remain bound to the membrane of the QIAamp MinElute column, contaminants are efficiently washed away using first Buffer AW1 and then Buffer AW2. For cleanup of genomic DNA, a single wash step with Buffer AW2 is sufficient to remove contaminants.

Eluting pure DNA

DNA is eluted from the QIAamp MinElute column using a small volume of either Buffer AE or distilled water. Since the elution volume is small, the eluted DNA will be concentrated.

QIAamp DNA Micro Procedure

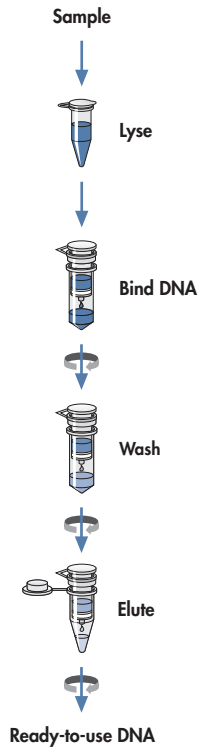

Automated purification

Purification of genomic DNA can be fully automated on the QIAcube®. The innovative QIAcube uses advanced technology to process QIAGEN spin columns, enabling seamless integration of automated, low-throughput sample prep into your laboratory workflow. Sample preparation using the QIAcube follows the same steps as the manual procedure (i.e., lyse, bind, wash, and elute), enabling you to continue using the QIAamp DNA Micro Kit for purification of high-quality DNA. For more information about the automated procedure, see the relevant protocol sheet available at www.qiagen.com/MyQIAcube.

The QIAcube is preinstalled with protocols for purification of plasmid DNA, genomic DNA, RNA, viral nucleic acids, and proteins, plus DNA and RNA cleanup. The range of protocols available is continually expanding, and additional QIAGEN protocols can be downloaded free of charge at www.qiagen.com/MyQIAcube.



The QIAcube.

Equipment and Reagents to Be Supplied by User

When working with chemicals, always wear a suitable lab coat, disposable gloves, and protective goggles. For more information, consult the appropriate material safety data sheets (MSDSs), available from the product supplier.

- Ethanol (96–100%)
- 1.5 ml or 2 ml microcentrifuge tubes (for lysis steps)
- 1.5 ml microcentrifuge tubes (for elution steps) (available from Brinkmann [Safe-Lock, cat. no. 022363204], Eppendorf [Safe-Lock, cat. no. 0030 120.086], or Sarstedt [Safety Cap, cat. no. 72.690])*
- Pipet tips (to avoid cross-contamination, we recommend pipet tips with aerosol barriers)
- Disposable gloves
- Thermomixer, heated orbital incubator, heating block, or water bath
- Microcentrifuge with rotor for 2 ml tubes
- Scissors
- Blood collection cards (if processing dried blood spots) (available from Whatman [903 paper or FTA® Card])*
- Paper punch for blood collection cards (available from Whatman)*
- Dithiothreitol (DTT) (if processing urine containing semen)

* This is not a complete list of suppliers and does not include many important vendors of biological supplies.

Important Notes

Eluting pure DNA

For downstream applications that require small starting volumes (e.g., some PCR assays), a more concentrated eluate may increase assay sensitivity. QIAamp MinElute columns allow a minimum elution volume of 20 μ l for concentrated nucleic acid eluates.

For downstream applications that require a larger starting volume, the elution volume can be increased to 100 μ l. However, an increase in elution volume will decrease the concentration of DNA in the eluate.

The volume of eluate recovered may be up to 5 μ l less than the volume of elution solution (i.e., Buffer AE or distilled water) applied to the QIAamp MinElute column. For example, an elution volume of 20 μ l results in \geq 15 μ l eluate. The volume of eluate recovered depends on the nature of the sample.

Buffer AE or distilled water should be equilibrated to room temperature (15–25°C) before it is applied to the QIAamp MinElute column. Yields will be increased if the column is incubated with Buffer AE or distilled water at room temperature for 5 minutes before centrifugation.

Eluted DNA can be collected in standard 1.5 ml microcentrifuge tubes (not provided). If the purified DNA is to be stored for up to 24 hours, we recommend storage at 2–8°C. For storage longer than 24 hours, we recommend elution in Buffer AE (10 mM Tris-Cl, 0.5 mM EDTA, pH 9.0) and storage at –20°C.

If high pH or EDTA affects sensitive downstream applications, use water for elution. However, ensure that the pH of the water is at least 7.0 (deionized water from certain sources can be acidic). DNA stored in water is subject to degradation by acid hydrolysis.

Yield and size of DNA

The yield of DNA isolated from biological samples strongly depends on the amount and source of the starting material. In addition, DNA size can also vary between 30 kb (in fresh samples) and 500 bp (in laser microdissections), or be below this range in forensic samples.

If DNA yields are below 1 μ g, quantification using a spectrophotometer will be difficult. In this case, we recommend quantitative amplification methods for determination of yield.

If purifying DNA from very small amounts of sample, such as 1–10 μ l blood or forensic samples, we recommend adding carrier RNA to the sample. Be aware that the sample will then contain considerably more carrier RNA than DNA. Therefore, when quantifying the purified DNA, avoid quantification methods that are not specific for DNA.

Carrier RNA

The kit is supplied with carrier RNA, which can be added to Buffer AL if required. Carrier RNA enhances binding of DNA to the QIAamp MinElute column membrane, especially if there are very few target molecules in the sample.

The amount of lyophilized carrier RNA provided is sufficient for the volume of Buffer AL supplied with the kit. The concentration of carrier RNA used in the QIAamp DNA Micro procedure allows the procedure to be used as a generic purification system compatible with many different amplification systems.

Different amplification systems vary in efficiency depending on the total amount of nucleic acid present in the reaction. If carrier RNA is used, eluates from QIAamp MinElute columns contain both sample DNA and carrier RNA, with the amount of carrier RNA greatly exceeding the amount of DNA. Calculations of how much eluate to add to downstream amplifications should therefore be based on the amount of carrier RNA added to Buffer AL. To obtain the highest levels of sensitivity in amplification reactions, it may be necessary to adjust the amount of carrier RNA added to Buffer AL.

Handling of QIAamp MinElute columns

Due to the sensitivity of nucleic acid amplification technologies, the following precautions are necessary when handling QIAamp MinElute columns, to avoid cross-contamination between sample preparations:

- Carefully apply the sample or solution to the QIAamp MinElute column. Pipet the sample into the QIAamp MinElute column without wetting the rim of the column.
- Always change pipet tips between liquid transfers. We recommend the use of aerosol-barrier pipet tips.
- Avoid touching the QIAamp MinElute column membrane with the pipet tip.
- After all pulse-vortexing steps, briefly centrifuge the microcentrifuge tubes to remove drops from the inside of the lids.
- Open only one QIAamp MinElute column at a time, and take care to avoid generating aerosols.
- Wear gloves throughout the entire procedure. In case of contact between gloves and sample, change gloves immediately.

Centrifugation

QIAamp MinElute columns will fit into most standard 1.5–2 ml microcentrifuge tubes. Additional 2 ml collection tubes are available separately.

Centrifugation of QIAamp MinElute columns is performed at 6000 × g (8000 rpm) to reduce centrifuge noise. Centrifugation at full speed will not improve DNA yields.

However, centrifugation of QIAamp MinElute columns at full speed is required in 2 steps of the procedure: in the dry centrifugation step after the membranes are washed, and in the elution step.

All centrifugation steps should be carried out at room temperature (15–25°C).

Processing QIAamp MinElute columns in a microcentrifuge

- Always close QIAamp MinElute columns before placing them in the microcentrifuge. Centrifuge as described in the relevant protocol.
- Flow-through fractions may contain hazardous waste and should be disposed of appropriately.
- For efficient parallel processing of multiple samples, we recommend filling a rack with collection tubes into which QIAamp MinElute columns can be transferred after centrifugation. Used collection tubes containing flow-through can be discarded, and the new collection tubes containing the QIAamp MinElute columns can be placed directly in the microcentrifuge.

Preparation of buffers

Preparing Buffer ATL

Before starting the procedure, check whether precipitate has formed in Buffer ATL. If necessary, dissolve by heating to 70°C with gentle agitation.

Preparing Buffer AL

Before starting the procedure, check whether precipitate has formed in Buffer AL. If necessary, dissolve by heating to 70°C with gentle agitation.

Preparing Buffer AW1

Add 25 ml ethanol (96–100%) to the bottle containing 19 ml Buffer AW1 concentrate. Tick the check box on the bottle label to indicate that ethanol has been added. Reconstituted Buffer AW1 can be stored at room temperature (15–25°C) for up to 1 year.

Note: Before starting the procedure, mix the reconstituted Buffer AW1 by shaking.

Preparing Buffer AW2

Add 30 ml ethanol (96–100%) to the bottle containing 13 ml Buffer AW2 concentrate. Reconstituted Buffer AW2 can be stored at room temperature (15–25°C) for up to 1 year.

Note: Before starting the procedure, mix the reconstituted Buffer AW2 by shaking.

Adding carrier RNA to Buffer AL

For purification of DNA from very small amounts of sample, such as low volumes of blood (<10 μ l), we recommend adding carrier RNA to Buffer AL. For samples containing larger amounts of DNA, addition of carrier RNA is optional.

Add 310 μ l Buffer AE to the tube containing 310 μ g lyophilized carrier RNA to obtain a solution of 1 μ g/ μ l. Dissolve the carrier RNA thoroughly, divide it into conveniently sized aliquots, and store at -20°C . Do not freeze–thaw the aliquots of carrier RNA more than 3 times.

Calculate the volume of Buffer AL and dissolved carrier RNA needed per batch of samples by multiplying the number of samples to be **simultaneously** processed by the volumes given in Table 1. To allow for pipetting errors, always prepare enough buffer for processing two extra samples.

Gently mix Buffer AL and dissolved carrier RNA by inverting the tube 10 times. To avoid foaming, do not vortex. Note that carrier RNA does not dissolve in Buffer AL. It must first be dissolved in Buffer AE and then added to Buffer AL. Buffer AL containing carrier RNA is stable at room temperature (15–25 $^{\circ}\text{C}$) for up to 48 hours.

Table 1. Volumes of Buffer AL and dissolved carrier RNA required for one DNA preparation using the QIAamp DNA Micro Kit

Protocol for isolation of genomic DNA from	Volume of Buffer AL added to sample	Dissolved carrier RNA
Small volumes of blood (page 16)	100 μ l	1 μ l
Dried blood spots (page 19)	200 μ l	1 μ l
Urine (page 22)	300 μ l	1 μ l
Tissues (page 25)	200 μ l	1 μ l
Laser-microdissected tissues (page 28)	50 μ l	1 μ l

Protocol: Isolation of Genomic DNA from Small Volumes of Blood

This protocol is for isolation of genomic DNA from 1–100 μl of whole blood treated with EDTA, citrate, or heparin-based anticoagulants.

Important point before starting

- Perform all centrifugation steps at room temperature (15–25°C).

Things to do before starting

- Equilibrate samples to room temperature (15–25°C).
- Equilibrate Buffer AE or distilled water for elution to room temperature.
- Set a thermomixer or heated orbital incubator to 56°C for use in step 5.
- Ensure that Buffers AW1 and AW2 have been prepared according to the instructions on page 14.
- If Buffer AL or Buffer ATL contains precipitates, dissolve by heating to 70°C with gentle agitation.
- Optional: If processing low volume samples (<10 μl), add carrier RNA dissolved in Buffer AE to Buffer AL, according to the instructions on page 15.

Procedure

1. **Pipet 1–100 μl whole blood into a 1.5 ml microcentrifuge tube (not provided).**
2. **Add Buffer ATL to a final volume of 100 μl .**
3. **Add 10 μl proteinase K.**
4. **Add 100 μl Buffer AL, close the lid, and mix by pulse-vortexing for 15 s.**

To ensure efficient lysis, it is essential that the sample, Buffer ATL, proteinase K, and Buffer AL are thoroughly mixed to yield a homogeneous solution.

Note: If the volume of blood is lower than 10 μl , we recommend adding carrier RNA to Buffer AL (see page 14). Note that carrier RNA does not dissolve in Buffer AL. It must first be dissolved in Buffer AE and then added to Buffer AL.

A white precipitate may form when Buffer AL is added to Buffer ATL. The precipitate does not interfere with the QIAamp procedure and will dissolve during the heat incubation in step 5.

5. **Incubate at 56°C for 10 min.**

Note: If samples are shaken during the incubation, DNA yields can be increased.

6. **Briefly centrifuge the 1.5 ml tube to remove drops from inside the lid.**

7. **Add 50 μ l ethanol (96–100%), close the lid, and mix thoroughly by pulse-vortexing for 15 s. Incubate for 3 min at room temperature.**

Note: If room temperature exceeds 25°C, cool the ethanol on ice before adding to the tube.

8. **Briefly centrifuge the 1.5 ml tube to remove drops from inside the lid.**
9. **Carefully transfer the entire lysate from step 8 to the QIAamp MinElute column (in a 2 ml collection tube) without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.**

If the lysate has not completely passed through the membrane after centrifugation, centrifuge again at a higher speed until the QIAamp MinElute column is empty.

10. **Carefully open the QIAamp MinElute column and add 500 μ l Buffer AW1 without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.**
11. **Carefully open the QIAamp MinElute column and add 500 μ l Buffer AW2 without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.**

Contact between the QIAamp MinElute column and the flow-through should be avoided. Some centrifuge rotors may vibrate upon deceleration, resulting in the flow-through, which contains ethanol, coming into contact with the QIAamp MinElute column. Take care when removing the QIAamp MinElute column and collection tube from the rotor, so that flow-through does not come into contact with the QIAamp MinElute column.

12. **Centrifuge at full speed (20,000 \times g; 14,000 rpm) for 3 min to dry the membrane completely.**

This step is necessary, since ethanol carryover into the eluate may interfere with some downstream applications.

13. **Place the QIAamp MinElute column in a clean 1.5 ml microcentrifuge tube (not provided) and discard the collection tube containing the flow-through. Carefully open the lid of the QIAamp MinElute column and apply 20–100 μ l Buffer AE or distilled water to the center of the membrane.**

If high pH or EDTA affects sensitive downstream applications, use water for elution (see page 12).

Important: Ensure that Buffer AE or distilled water is equilibrated to room temperature (15–25°C). If using small elution volumes (<50 μ l), dispense Buffer AE or distilled water onto the center of the membrane to ensure complete elution of bound DNA.

QIAamp MinElute columns provide flexibility in the choice of elution volume. Choose a volume according to the requirements of the downstream application. Remember that the volume of eluate will be up to 5 μl less than the volume of elution solution applied to the column.

- 14. Close the lid and incubate at room temperature (15–25°C) for 1 min. Centrifuge at full speed (20,000 \times g; 14,000 rpm) for 1 min.**

Incubating the QIAamp MinElute column loaded with Buffer AE or water for 5 min at room temperature before centrifugation generally increases DNA yield.

Protocol: Isolation of Genomic DNA from Dried Blood Spots

This protocol is for isolation of genomic DNA from blood card punches with diameters of approximately 3 mm (1/8 inch). This protocol is suitable for untreated blood or blood treated with anticoagulants such as EDTA, citrate, or heparin. The blood must be spotted and dried on filter paper such as 903 paper or FTA Card (Whatman), Guthrie test cards, or comparable blood cards.

Important point before starting

- Perform all centrifugation steps at room temperature (15–25°C).

Things to do before starting

- Equilibrate Buffer AE or distilled water for elution to room temperature (15–25°C).
- If Buffer AL or Buffer ATL contains precipitates, dissolve by heating to 70°C with gentle agitation.
- Set a thermomixer or heated orbital incubator to 56°C for use in step 4, and a second thermomixer or heated orbital incubator to 70°C for use in step 7. If thermomixers or heated orbital incubators are not available, heating blocks or water baths can be used instead.
- Ensure that Buffers AW1 and AW2 have been prepared according to the instructions on page 14.
- Optional: If processing very small amounts of starting material, add carrier RNA dissolved in Buffer AE to Buffer AL, according to the instructions on page 15.

Procedure

1. **Cut 3 mm (1/8 inch) diameter punches from a dried blood spot with a single-hole paper punch. Place up to 3 blood card punches into a 1.5 ml microcentrifuge tube (not provided).**
2. **Add 180 µl Buffer ATL.**
3. **Add 20 µl proteinase K and mix thoroughly by vortexing.**
4. **Place the 1.5 ml tube in a thermomixer or heated orbital incubator, and incubate at 56°C with shaking at 900 rpm for 1 h.**
If using a heating block or water bath, vortex the tube for 10 s every 10 min to improve lysis.
5. **Briefly centrifuge the 1.5 ml tube to remove drops from inside the lid.**

6. Add 200 μ l Buffer AL, close the lid, and mix by pulse-vortexing for 10 s.

To ensure efficient lysis, it is essential that the sample and Buffer AL are thoroughly mixed to yield a homogeneous solution.

A white precipitate may form when Buffer AL is added to Buffer ATL. The precipitate does not interfere with the QIAamp procedure and will dissolve during the heat incubation in step 7.

Note: If processing only 1 blood card punch with a diameter of 3 mm or less, we recommend adding carrier RNA to Buffer AL (see page 14). Note that carrier RNA does not dissolve in Buffer AL. It must first be dissolved in Buffer AE and then added to Buffer AL.

7. Place the 1.5 ml tube in a thermomixer or heated orbital incubator, and incubate at 70°C for 10 min, with shaking at 900 rpm.

If using a heating block or water bath, vortex the tube for 10 s every 3 min.

8. Briefly centrifuge the 1.5 ml tube to remove drops from inside the lid.**9. Carefully transfer the entire lysate from step 8 to the QIAamp MinElute column (in a 2 ml collection tube) without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.**

If the lysate has not completely passed through the membrane after centrifugation, centrifuge again at a higher speed until the QIAamp MinElute column is empty.

10. Carefully open the QIAamp MinElute column and add 500 μ l Buffer AW1 without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.**11. Carefully open the QIAamp MinElute column and add 500 μ l Buffer AW2 without wetting the rim. Close the lid and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.**

Contact between the QIAamp MinElute column and the flow-through should be avoided. Some centrifuge rotors may vibrate upon deceleration, resulting in the flow-through, which contains ethanol, coming into contact with the QIAamp MinElute column. Take care when removing the QIAamp MinElute column and collection tube from the rotor, so that flow-through does not come into contact with the QIAamp MinElute column.

12. Centrifuge at full speed (20,000 \times g; 14,000 rpm) for 3 min to dry the membrane completely.

This step is necessary, since ethanol carryover into the eluate may interfere with some downstream applications.

13. Place the QIAamp MinElute column in a clean 1.5 ml microcentrifuge tube (not provided) and discard the collection tube containing the flow-through. Carefully open the lid of the QIAamp MinElute column and apply 20–100 μ l Buffer AE or distilled water to the center of the membrane.

If high pH or EDTA affects sensitive downstream applications, use water for elution (see page 12).

Important: Ensure that Buffer AE or distilled water is equilibrated to room temperature (15–25°C). If using small elution volumes (<50 μ l), dispense Buffer AE or distilled water onto the center of the membrane to ensure complete elution of bound DNA.

QIAamp MinElute columns provide flexibility in the choice of elution volume. Choose a volume according to the requirements of the downstream application. Remember that the volume of eluate will be up to 5 μ l less than the volume of elution solution applied to the column.

14. Close the lid and incubate at room temperature (15–25°C) for 1 min. Centrifuge at full speed (20,000 \times g; 14,000 rpm) for 1 min.

Incubating the QIAamp MinElute column loaded with Buffer AE or water for 5 min at room temperature before centrifugation generally increases DNA yield.

Protocol: Isolation of Genomic DNA from Urine

This protocol is for isolation of genomic DNA from urine.

Important points before starting

- Perform all centrifugation steps at room temperature (15–25°C).
- Check whether carrier RNA is required (see pages 13 and 15).

Things to do before starting

- Equilibrate samples to room temperature (15–25°C).
- Equilibrate Buffer AE or distilled water for elution to room temperature.
- Set a thermomixer or heated orbital incubator to 56°C for use in step 5. If a thermomixer or heated orbital incubator is not available, a heating block or water bath can be used instead.
- If Buffer AL or Buffer ATL contains precipitates, dissolve by heating to 70°C with gentle agitation.
- Optional: Prepare an aqueous 1 M DTT (dithiothreitol) stock solution. Store aliquots at –20°C. Thaw immediately before use.
- Ensure that Buffers AW1 and AW2 have been prepared according to the instructions on page 14.

Procedure

1. **Transfer up to 1 ml urine to a 1.5 ml microcentrifuge tube (not provided) and centrifuge at 6000 x g (8000 rpm) for 2 min.**

If the volume of the urine sample is between 1 and 10 ml, pellet the cells in an appropriately sized centrifugation tube. Discard the supernatant, add 500 µl Buffer AE to the pellet, and vortex for 5 s. Briefly centrifuge and transfer the sample to a 1.5 ml microcentrifuge tube (not provided). Continue the procedure from step 3.

2. **Discard the supernatant, add 500 µl Buffer AE, and vortex for 5 s.**
3. **Centrifuge at 6000 x g (8000 rpm) for 2 min.**
4. **Discard the supernatant, add 300 µl Buffer ATL and 20 µl proteinase K to the pellet, and mix by pulse-vortexing for 10 s.**

Adding 20 µl 1 M DTT may increase sensitivity, since urine can contain sperm cells which can only be lysed in the presence of reducing agents such as DTT or β-mercaptoethanol.

5. **Place the 1.5 ml tube in a thermomixer or heated orbital incubator, and incubate at 56°C, with shaking at 900 rpm, for 1 h.**

If using a heating block or water bath, vortex the tube for 10 s every 15 min to improve lysis.

6. **Briefly centrifuge the 1.5 ml tube to remove drops from inside the lid.**
7. **Add 300 μ l Buffer AL and 50 μ l ethanol (96–100%), close the lid, and mix by pulse-vortexing for 10 s.**

To ensure efficient binding in step 9, it is essential that the sample, Buffer AL, and ethanol are thoroughly mixed to yield a homogeneous solution.

Note: We recommend adding carrier RNA to Buffer AL (see page 14). Note that carrier RNA does not dissolve in Buffer AL. It must first be dissolved in Buffer AE and then added to Buffer AL.

8. **Briefly centrifuge the 1.5 ml tube to remove drops from inside the lid.**

Note: The white precipitate which may appear in step 7 does not need to be pelleted. It can be transferred to the QIAamp MinElute column together with the lysate, and does not interfere with the QIAamp procedure.

9. **Carefully transfer the supernatant from step 8 to the QIAamp MinElute column (in a 2 ml collection tube) without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.**

If the lysate has not completely passed through the membrane after centrifugation, centrifuge again at a higher speed until the QIAamp MinElute column is empty.

10. **Carefully open the QIAamp MinElute column and add 500 μ l Buffer AW1 without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.**
11. **Carefully open the QIAamp MinElute column and add 500 μ l Buffer AW2 without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.**

Contact between the QIAamp MinElute column and the flow-through should be avoided. Some centrifuge rotors may vibrate upon deceleration, resulting in the flow-through, which contains ethanol, coming into contact with the QIAamp MinElute column. Take care when removing the QIAamp MinElute column and collection tube from the rotor, so that flow-through does not come into contact with the QIAamp MinElute column.

- 12. Centrifuge at full speed (20,000 x g; 14,000 rpm) for 3 min to dry the membrane completely.**

This step is necessary, since ethanol carryover into the eluate may interfere with some downstream applications.

- 13. Place the QIAamp MinElute column in a clean 1.5 ml microcentrifuge tube (not provided) and discard the collection tube containing the flow-through. Carefully open the lid of the QIAamp MinElute column and apply 20–50 μ l Buffer AE or distilled water to the center of the membrane.**

If high pH or EDTA affects sensitive downstream applications, use water for elution (see page 12).

Important: Ensure that Buffer AE or distilled water is equilibrated to room temperature (15–25°C). Dispense Buffer AE or distilled water onto the center of the membrane to ensure complete elution of bound DNA.

QIAamp MinElute columns provide flexibility in the choice of elution volume. Choose a volume according to the requirements of the downstream application. Remember that the volume of eluate will be up to 5 μ l less than the volume of elution solution applied to the column.

- 14. Close the lid and incubate at room temperature (15–25°C) for 1 min. Centrifuge at full speed (20,000 x g; 14,000 rpm) for 1 min.**

Incubating the QIAamp MinElute column loaded with Buffer AE or water for 5 min at room temperature before centrifugation generally increases DNA yield.

Protocol: Isolation of Genomic DNA from Tissues

This protocol is for isolation of genomic DNA from less than 10 mg tissue.

Important points before starting

- Perform all centrifugation steps at room temperature (15–25°C)
- If isolating DNA from very small amounts of tissue, carrier RNA is required (see page 15).
- Prepare tissue samples on a cold surface (e.g., a glass, steel, or aluminum plate placed on top of a block of dry ice).
- If using frozen tissue, ensure that the sample does not thaw out before addition of Buffer ATL in step 2.

Things to do before starting

- Equilibrate Buffer AE or distilled water for elution to room temperature (15–25°C).
- Set a thermomixer or heated orbital incubator to 56°C for use in step 4. If a thermomixer or heated orbital incubator is not available, a heating block or water bath can be used instead.
- If Buffer AL or Buffer ATL contains precipitates, dissolve by heating to 70°C with gentle agitation.
- Ensure that Buffers AW1 and AW2 have been prepared according to the instructions on page 14.

Procedure

1. **Transfer a tissue sample of less than 10 mg in weight to a 1.5 ml microcentrifuge tube (not provided).**
2. **Immediately add 180 μ l Buffer ATL, and equilibrate to room temperature (15–25°C).**
3. **Add 20 μ l proteinase K and mix by pulse-vortexing for 15 s.**
4. **Place the 1.5 ml tube in a thermomixer or heated orbital incubator, and incubate at 56°C overnight until the sample is completely lysed.**

For small amounts of tissue, lysis is complete in 4–6 h, but best results are achieved after overnight lysis.

5. Add 200 μ l Buffer AL, close the lid, and mix by pulse-vortexing for 15 s.

To ensure efficient lysis, it is essential that the sample and Buffer AL are thoroughly mixed to yield a homogenous solution.

Note: If carrier RNA is required (see page 13), add 1 μ g dissolved carrier RNA to 200 μ l Buffer AL. Note that carrier RNA does not dissolve in Buffer AL. It must first be dissolved in Buffer AE and then added to Buffer AL.

6. Add 200 μ l ethanol (96–100%), close the lid, and mix thoroughly by pulse-vortexing for 15 s. Incubate for 5 min at room temperature (15–25°C).

Note: If room temperature exceeds 25°C, cool the ethanol on ice before adding to the tube.

7. Briefly centrifuge the 1.5 ml tube to remove drops from inside the lid.

8. Carefully transfer the entire lysate from step 7 to the QIAamp MinElute column (in a 2 ml collection tube) without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.

If the lysate has not completely passed through the membrane after centrifugation, centrifuge again at a higher speed until the QIAamp MinElute column is empty.

9. Carefully open the QIAamp MinElute column and add 500 μ l Buffer AW1 without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.

10. Carefully open the QIAamp MinElute column and add 500 μ l Buffer AW2 without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.

Contact between the QIAamp MinElute column and the flow-through should be avoided. Some centrifuge rotors may vibrate upon deceleration, resulting in the flow-through, which contains ethanol, coming into contact with the QIAamp MinElute column. Take care when removing the QIAamp MinElute column and collection tube from the rotor, so that flow-through does not come into contact with the QIAamp MinElute column.

11. Centrifuge at full speed (20,000 \times g; 14,000 rpm) for 3 min to dry the membrane completely.

This step is necessary, since ethanol carryover into the eluate may interfere with some downstream applications.

- 12. Place the QIAamp MinElute column in a clean 1.5 ml microcentrifuge tube (not provided) and discard the collection tube containing the flow-through. Carefully open the lid of the QIAamp MinElute column and apply 20–100 μ l Buffer AE or distilled water to the center of the membrane.**

If high pH or EDTA affects sensitive downstream applications, use water for elution (see page 12).

Important: Ensure that Buffer AE or distilled water is equilibrated to room temperature (15–25°C). If using small elution volumes (<50 μ l), dispense Buffer AE or distilled water onto the center of the membrane to ensure complete elution of bound DNA.

QIAamp MinElute columns provide flexibility in the choice of elution volume. Choose a volume according to the requirements of the downstream application. Remember that the volume of eluate will be up to 5 μ l less than the volume of the solution applied to the column.

- 13. Close the lid and incubate at room temperature (15–25°C) for 1 min. Centrifuge at full speed (20,000 \times g; 14,000 rpm) for 1 min.**

Incubating the QIAamp MinElute column loaded with Buffer AE or water for 5 min at room temperature before centrifugation generally increases DNA yield.

Protocol: Isolation of Genomic DNA from Laser-Microdissected Tissues

This protocol is for isolation of genomic DNA from laser-microdissected tissue. Laser-microdissected tissue specimens present a particular challenge for molecular analysis, as nucleic acids must be purified from very small amounts of starting material. In addition, fixation and staining steps may compromise the integrity of DNA, and it may be necessary either to modify fixation protocols or to use cryosections from flash-frozen specimens to minimize this problem.

A wide range of equipment and consumables for sectioning, staining, and microdissection of specimens is available from Leica (www.leica-microsystems.com).

Important points before starting

- Perform all centrifugation steps at room temperature (15–25°C)
- If isolating DNA from very small numbers of cells, carrier RNA is required (see page 15).

Things to do before starting

- Equilibrate samples to room temperature (15–25°C).
- Equilibrate Buffer AE or distilled water for elution to room temperature.
- Set a thermomixer or heated orbital incubator to 56°C for use in step 3. If a thermomixer or heated orbital incubator is not available, a heating block or water bath can be used instead.
- If Buffer AL or Buffer ATL contains precipitates, dissolve by heating to 70°C with gentle agitation.
- Ensure that Buffers AW1 and AW2 have been prepared according to the instructions on page 14.

Procedure

1. **Add 15 µl Buffer ATL to a laser-microdissected sample collected in a 0.2 ml microcentrifuge tube (not provided).**
2. **Add 10 µl proteinase K and mix by pulse-vortexing for 15 s.**
3. **Place the 0.2 ml tube in a thermomixer or heated orbital incubator, and incubate at 56°C for 3 h (16 h for formalin-fixed tissues) with occasional agitation.**
The incubation time may vary depending on the amount of tissue collected.
4. **Add 25 µl Buffer ATL.**

5. Add 50 μ l Buffer AL, close the lid, and mix by pulse-vortexing for 15 s.

To ensure efficient lysis, it is essential that the sample and Buffer AL are thoroughly mixed to yield a homogenous solution.

6. Add 50 μ l ethanol (96–100%), close the lid, and mix thoroughly by pulse-vortexing for 15 s. Incubate for 5 min at room temperature (15–25°C).

Note: If room temperature exceeds 25°C, cool the ethanol on ice before adding to the tube.

7. Briefly centrifuge the 0.2 ml tube to remove drops from inside the lid.

8. Carefully transfer the entire lysate from step 7 to the QIAamp MinElute column (in a 2 ml collection tube) without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.

If the lysate has not completely passed through the membrane after centrifugation, centrifuge again at a higher speed until the QIAamp MinElute column is empty.

9. Carefully open the QIAamp MinElute column and add 500 μ l Buffer AW1 without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.

10. Carefully open the QIAamp MinElute column and add 500 μ l Buffer AW2 without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.

Contact between the QIAamp MinElute column and the flow-through should be avoided. Some centrifuge rotors may vibrate upon deceleration, resulting in the flow-through, which contains ethanol, coming into contact with the QIAamp MinElute column. Take care when removing the QIAamp MinElute column and collection tube from the rotor, so that flow-through does not come into contact with the QIAamp MinElute column.

11. Centrifuge at full speed (20,000 \times g; 14,000 rpm) for 3 min to dry the membrane completely.

This step is necessary, since ethanol carryover into the eluate may interfere with some downstream applications.

12. Place the QIAamp MinElute column in a clean 1.5 ml microcentrifuge tube (not provided) and discard the collection tube containing the flow-through. Carefully open the lid of the QIAamp MinElute column and apply 20–100 μ l Buffer AE or distilled water to the center of the membrane.

If high pH or EDTA affects sensitive downstream applications, use water for elution (see page 12).

Important: Ensure that Buffer AE or distilled water is equilibrated to room temperature (15–25°C). Dispense Buffer AE or distilled water onto the center of the membrane to ensure complete elution of bound DNA.

QIAamp MinElute columns provide flexibility in the choice of elution volume. Choose a volume according to the requirements of the downstream application. Remember that the volume of eluate will be up to 5 μ l less than the volume of the solution applied to the column.

- 13. Close the lid and incubate at room temperature (15–25°C) for 1 min. Centrifuge at full speed (20,000 x g; 14,000 rpm) for 1 min.**

Incubating the QIAamp MinElute column loaded with Buffer AE or water for 5 min at room temperature before centrifugation generally increases DNA yield.

Protocol: Cleanup of Genomic DNA

This protocol is for cleanup of genomic DNA. Use this protocol to restore the suitability of the DNA for PCR, or to increase the concentration of the DNA.

Important point before starting

- Perform all centrifugation steps at room temperature (15–25°C).

Things to do before starting

- Equilibrate samples to room temperature (15–25°C).
- Equilibrate Buffer AE or distilled water for elution to room temperature.
- Ensure that Buffers AW1 and AW2 have been prepared according to the instructions on page 14

Procedure

- 1. Add up to 100 μ l of genomic DNA (containing up to 10 μ g DNA) to a 1.5 ml microcentrifuge tube (not provided).**
If the volume of DNA is less than 100 μ l, add deionized water to a final volume of 100 μ l.
- 2. Add 10 μ l Buffer AW1.**
- 3. Add 250 μ l Buffer AW2 and mix by pulse-vortexing for 10 s.**
- 4. Transfer the entire sample from step 3 to the QIAamp MinElute column without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.**
- 5. Carefully open the QIAamp MinElute column and add 500 μ l Buffer AW2 without wetting the rim. Close the lid, and centrifuge at 6000 \times g (8000 rpm) for 1 min. Place the QIAamp MinElute column in a clean 2 ml collection tube, and discard the collection tube containing the flow-through.**

Contact between the QIAamp MinElute column and the flow-through should be avoided. Some centrifuge rotors may vibrate upon deceleration, resulting in the flow-through, which contains ethanol, coming into contact with the QIAamp MinElute column. Take care when removing the QIAamp MinElute column and collection tube from the rotor, so that flow-through does not come into contact with the QIAamp MinElute column.

- 6. Centrifuge at full speed (20,000 x g; 14,000 rpm) for 3 min to dry the membrane completely.**

This step is necessary, since ethanol carryover into the eluate may interfere with some downstream applications.

- 7. Place the QIAamp MinElute column in a clean 1.5 ml microcentrifuge tube (not provided) and discard the collection tube containing the flow-through. Carefully open the lid of the QIAamp MinElute column and apply 20–100 μ l Buffer AE or distilled water to the center of the membrane.**

If high pH or EDTA affects sensitive downstream applications, use water for elution (see page 12).

Important: Ensure that Buffer AE or distilled water is equilibrated to room temperature (15–25°C). If using small elution volumes (<50 μ l), dispense Buffer AE or distilled water onto the center of the membrane to ensure complete elution of bound DNA.

QIAamp MinElute columns provide flexibility in the choice of elution volume. Choose a volume according to the requirements of the downstream application. Remember that the volume of eluate will be up to 5 μ l less than the volume of elution solution applied to the column.

- 8. Close the lid and incubate at room temperature (15–25°C) for 1 min. Centrifuge at full speed (20,000 x g; 14,000 rpm) for 1 min.**

Incubating the QIAamp MinElute column loaded with Buffer AE or water for 5 min at room temperature before centrifugation generally increases DNA yield.

Troubleshooting Guide

This troubleshooting guide may be helpful in solving any problems that may arise. For more information, see also the Frequently Asked Questions page at our Technical Support Center: www.qiagen.com/FAQ/FAQlist.aspx. The scientists in QIAGEN Technical Services are always happy to answer any questions you may have about either the information and protocols in this handbook or sample and assay technologies (for contact information, see back cover or visit www.qiagen.com).

Comments and suggestions

Little or no DNA in the eluate

- | | |
|---|---|
| a) Carrier RNA was not added to Buffer AL | Dissolve carrier RNA in Buffer AE and mix with Buffer AL as described on page 15. Repeat the purification procedure with new samples. |
| b) Samples were frozen and thawed more than once | Avoid repeated freezing and thawing of samples. Where possible, always use fresh samples or samples that have been thawed only once. |
| c) Samples were kept at room temperature for too long | DNA in the samples may degrade during prolonged storage at room temperature. Where possible, always use fresh samples, or store the samples at 2–8°C (non-dried blood) or at –20°C (tissue samples). Dried blood spots or stains can be stored at room temperature in the dark without significant DNA degradation. |
| d) Insufficient sample lysis in Buffer AL | Proteinase K was stored at high temperatures for a prolonged time. Repeat the procedure using new samples and fresh proteinase K. |
| e) Buffer AL–carrier RNA mixture was mixed insufficiently | Mix Buffer AL with carrier RNA by gently inverting the tube at least 10 times. |
| f) Low-percentage ethanol was used instead of 96–100% | Repeat the purification procedure with new samples and 96–100% ethanol. |
| g) Buffer AW1 or AW2 was prepared incorrectly | Check that the Buffer AW1 and Buffer AW2 concentrates were diluted with the correct volume of 96–100% ethanol. Repeat the purification procedure with new samples, if available. |
| h) pH of water used for elution was too low | DNA does not dissolve easily in acidic solutions. Ensure that the pH of the water used for elution is >7.0. |

Comments and suggestions

DNA does not perform well in downstream enzymatic reactions

- | | | |
|----|--|---|
| a) | Little or no DNA in the eluate | See "Little or no DNA in the eluate" (page 43) for possible reasons. Increase the amount of eluate added to the reaction, if possible. |
| b) | Too much or too little carrier RNA in the eluate | Determine the maximum amount of carrier RNA suitable for your amplification reaction. Adjust the concentration of carrier RNA added to Buffer AL accordingly. |
| c) | Reduced sensitivity | Determine the maximum volume of eluate suitable for your amplification reaction. Reduce or increase the volume of eluate added to the amplification reaction accordingly. The elution volume can be adjusted proportionally.

Use water instead of Buffer AE for elution (see page 12). |
| d) | Performance of purified DNA in downstream assays varies with the age of the reconstituted wash buffers | Salt and ethanol components of wash Buffers AW1 and AW2 may have separated out after being unused for a long period. Always mix buffers thoroughly before each purification procedure. |
| e) | Eluate used for PCR contains too much Buffer AE | Buffer AE contains EDTA and is buffered at pH 9.0. Some <i>Taq</i> DNA Polymerases are inhibited by EDTA and/or have a narrow pH optimum. Use distilled water for elution instead of Buffer AE.

Use another DNA polymerase, such as QIAGEN HotStarTaq® DNA Polymerase. |

General handling

- | | | |
|----|--------------------------------|---|
| a) | Clogged QIAamp MinElute column | Incomplete lysis caused clogging of the membrane. Increase the lysis time to fully lyse the sample. |
| b) | Variable elution volumes | Different sample types have been processed. |

Appendix: Working with DNA

General handling

Proper microbiological aseptic technique should always be used when working with small sample sizes. Hands and dust particles may carry bacteria and molds, and are the most common sources of contamination. Always wear latex or vinyl gloves while handling reagents and samples to prevent contamination from the surface of the skin or from dusty laboratory equipment. Change gloves frequently and keep tubes closed.

Disposable plasticware

The use of sterile, disposable polypropylene tubes is recommended throughout the purification procedure. These tubes are generally DNase-free.

References

QIAGEN maintains a large, up-to-date online database of scientific publications utilizing QIAGEN products. Comprehensive search options allow you to find the articles you need, either by a simple keyword search or by specifying the application, research area, title, etc.

For a complete list of references, visit the QIAGEN Reference Database online at www.qiagen.com/RefDB/search.asp or contact QIAGEN Technical Services or your local distributor.

Ordering Information

Product	Contents	Cat. no.
QIAamp DNA Micro Kit (50)	For 50 DNA preps: 50 QIAamp MinElute Columns, Proteinase K, Carrier RNA, Buffers, Collection Tubes (2 ml)	56304
Accessories		
Buffer ATL (200 ml)	200 ml Tissue Lysis Buffer	19076
Buffer AL (216 ml)	216 ml Lysis Buffer	19075
Buffer AW1 (concentrate, 242 ml)	242 ml Wash Buffer (1) Concentrate	19081
Buffer AW2 (concentrate, 324 ml)	324 ml Wash Buffer (2) Concentrate	19072
Buffer AE (240 ml)	240 ml Elution Buffer	19077
QIAGEN Proteinase K (2 ml)*	2 ml (>600 mAU/ml, solution)	19131
Collection Tubes (2 ml)	1000 Collection Tubes (2 ml)	19201
Related products		
QIAamp DNA Blood Mini Kit — for isolation of DNA from blood and related body fluids		
QIAamp DNA Blood Mini Kit (50)*	For 50 DNA minipreps: 50 QIAamp Mini Spin Columns, QIAGEN Protease, Reagents, Buffers, Collection Tubes (2 ml)	51104
QIAamp DNA Mini Kit — for isolation of DNA from tissue, blood, and body fluids		
QIAamp DNA Mini Kit (50)*	For 50 DNA preps: 50 QIAamp Mini Spin Columns, Proteinase K, Reagents, Buffers, Collection Tubes (2 ml)	51304

For up-to-date licensing information and product-specific disclaimers, see the respective QIAGEN kit handbook or user manual. QIAGEN kit handbooks and user manuals are available at www.qiagen.com or can be requested from QIAGEN Technical Services or your local distributor.

* Larger kit sizes available; please inquire.

Notes

Notes

Trademarks: QIAGEN®, QIAamp®, QIAcube®, HotStarTaq®, MinElute™ (QIAGEN Group); FTA® and Whatman® (Whatman PLC).

Limited License Agreement

Use of this product signifies the agreement of any purchaser or user of the QIAamp DNA Micro Kit to the following terms:

1. The QIAamp DNA Micro Kit may be used solely in accordance with the *QIAamp DNA Micro Handbook* and for use with components contained in the Kit only. QIAGEN grants no license under any of its intellectual property to use or incorporate the enclosed components of this Kit with any components not included within this Kit except as described in the *QIAamp DNA Micro Handbook* and additional protocols available at www.qiagen.com.
2. Other than expressly stated licenses, QIAGEN makes no warranty that this Kit and/or its use(s) do not infringe the rights of third-parties.
3. This Kit and its components are licensed for one-time use and may not be reused, refurbished, or resold.
4. QIAGEN specifically disclaims any other licenses, expressed or implied other than those expressly stated.
5. The purchaser and user of the Kit agree not to take or permit anyone else to take any steps that could lead to or facilitate any acts prohibited above. QIAGEN may enforce the prohibitions of this Limited License Agreement in any Court, and shall recover all its investigative and Court costs, including attorney fees, in any action to enforce this Limited License Agreement or any of its intellectual property rights relating to the Kit and/or its components.

For updated license terms, see www.qiagen.com.

© 2010 QIAGEN, all rights reserved.

www.qiagen.com

Australia = Orders 1-800-243-800 = Fax 03-9840-9888 = Technical 1-800-243-066

Austria = Orders 0800-28-10-10 = Fax 0800/28-10-19 = Technical 0800-28-10-11

Belgium = Orders 0800-79612 = Fax 0800-79611 = Technical 0800-79556

Brazil = Orders 0800-557779 = Fax 55-11-5079-4001 = Technical 0800-557779

Canada = Orders 800-572-9613 = Fax 800-713-5951 = Technical 800-DNA-PREP (800-362-7737)

China = Orders 86-21-3865-3865 = Fax 86-21-3865-3965 = Technical 800-988-0325

Denmark = Orders 80-885945 = Fax 80-885944 = Technical 80-885942

Finland = Orders 0800-914416 = Fax 0800-914415 = Technical 0800-914413

France = Orders 01-60-920-926 = Fax 01-60-920-925 = Technical 01-60-920-930 = Offers 01-60-920-928

Germany = Orders 02103-29-12000 = Fax 02103-29-22000 = Technical 02103-29-12400

Hong Kong = Orders 800 933 965 = Fax 800 930 439 = Technical 800 930 425

Ireland = Orders 1800 555 049 = Fax 1800 555 048 = Technical 1800 555 061

Italy = Orders 800-789-544 = Fax 02-334304-826 = Technical 800-787980

Japan = Telephone 03-6890-7300 = Fax 03-5547-0818 = Technical 03-6890-7300

Korea (South) = Orders 080-000-7146 = Fax 02-2626-5703 = Technical 080-000-7145

Luxembourg = Orders 8002-2076 = Fax 8002-2073 = Technical 8002-2067

Mexico = Orders 01-800-7742-639 = Fax 01-800-1122-330 = Technical 01-800-7742-639

The Netherlands = Orders 0800-0229592 = Fax 0800-0229593 = Technical 0800-0229602

Norway = Orders 800-18859 = Fax 800-18817 = Technical 800-18712

Singapore = Orders 1800-742-4362 = Fax 65-6854-8184 = Technical 1800-742-4368

Spain = Orders 91-630-7050 = Fax 91-630-5145 = Technical 91-630-7050

Sweden = Orders 020-790282 = Fax 020-790582 = Technical 020-798328

Switzerland = Orders 055-254-22-11 = Fax 055-254-22-13 = Technical 055-254-22-12

UK = Orders 01293-422-911 = Fax 01293-422-922 = Technical 01293-422-999

USA = Orders 800-426-8157 = Fax 800-718-2056 = Technical 800-DNA-PREP (800-362-7737)

